

Gradian Health Systems®



Universal Anaesthesia Machine®
Maintenance Manual

HOW TO GET HELP	1
GENERAL DESCRIPTION OF THE UNIVERSAL ANAESTHESIA MACHINE®	1-2
SPECIFICATIONS	2-5
COMPONENTS OF THE UAM®	6-10
UNPACKING AND PREPARING FOR USE	11-15
RECOMMENDED UAM MAINTENANCE SCHEDULE	15
FULL FUNCTION TEST	16-17
MAINTENANCE PROCEDURES	17
Cleaning or Replacing the Oxygen Concentrator Filter	17-18
Replacing the Oxygen Sensor	18
Removing and Replacing the Draw-Over Vaporizer	18-19
Cleaning the Halothane Vaporizer	19
Replacing the Fenton Balloon	20
Replacing the Control Screen Circuit Board	21
Replacing the Concentrator Loss of Power Alarm Battery	22
Removing the Concentrator	23-24
Installing the Concentrator	25
Replacing the Control Screen Battery	26
Removing, Inspecting and Cleaning the Pressure Relief Valve	27-30
PROBLEM SOLVING	31
MAIN PARTS	32-33
ILLUSTRATED PARTS	33-37
APPENDIX I: BREATHING CIRCUIT DIAGRAM	38
APPENDIX II: OXYGEN CONCENTRATOR FLOW DIAGRAM	39
APPENDIX III: ELECTRICAL DIAGRAM	40
APPENDIX IV: UAM PREVENTIVE MAINTENANCE AND FUNCTIONAL TESTS FORM	41
APPENDIX V: UAM INSPECTION LABEL	42
APPENDIX VI: PRESSURE CONVERSIONS & COLOUR IDENTIFICATION FOR MEDICAL GASES	42-43

Copyright © 2016 by Gradian Health Systems, LLC. All rights reserved.

Gradian Health Systems, UAM and the Gradian logo are registered trademarks (*) of Gradian Health Systems, LLC, registered in the United States.

Printed in the United States of America. Rev 2.0, P/N: 1974-005, CE 0120.

The UAM is manufactured by OES Medical Ltd., Witney, Oxon, UK.



HOW TO GET HELP

EMAIL: service@gradianhealth.org
WEB: www.gradianhealth.org
PHONE: +1.212.537.0340
WHATSAPP: +1.929.280.0210

ADDRESS: Gradian Health Systems, 915 Broadway, Suite 1001, New York, NY 10010, USA

NOTE: For the latest UAM specifications and technical materials please visit us online at:
<http://www.gradianhealth.org/universal-anaesthesia-machine/>

DESCRIPTION OF THE UNIVERSAL ANAESTHESIA MACHINE®

The UAM® combines elements from both continuous flow and draw-over designs to provide a workstation that delivers controlled anaesthesia in all hospital settings.

The breathing system is a low pressure (limited to 5 cm H₂O) continuous flow system that defaults to a demand-flow (draw-over) system when the supply fails. There are three one-way valves and a balloon operated inflating valve which provide this uniquely versatile method of anaesthesia delivery.

Oxygen can be provided from a variety of sources, and if none of those sources is available the machine will automatically draw in room air without the need for the user to change settings. The vaporizer provides calibrated delivery of anaesthetic agent and the manual bellows ensures that safe assisted respiration can continue without requiring a high pressure gas source.

The system has been designed to accept standard anaesthesia accessories such as the breathing circuit, oxygen sensor, and standard connectors for the oxygen and nitrous oxide supplies. The UAM is also designed for simple field service. Most components can be removed with standard tools and local technicians can install replacement parts without removing the workstation from the hospital.

The UAM carries the CE Mark, having passed EU regulatory inspections. It is manufactured in an ISO-certified factory in England. It is used in a wide variety of hospitals, from internationally recognized academic centres to resource-constrained district hospitals.

The Universal Anaesthesia Machine ensures that an adequate supply of oxygen and anaesthesia is always available to meet patient demand, no matter what the environment.

The UAM provides several ways to deliver oxygen, including connections for cylinder and pipeline. Most importantly, the machine produces its own oxygen using an integrated concentrator that converts room air into 95% oxygen. Air is drawn into the UAM and is purified by 2 different filters. Ambient air is made up of 78% nitrogen, 21% oxygen and 1% other gases. The electrically powered oxygen concentrator in the UAM removes nitrogen from room air, creating a mixture that is 95% oxygen and 5% other gases.

The compressed air passes into a canister filled with a powdered material called Zeolite. Nitrogen sticks to the Zeolite, while oxygen and other gases flow through it to a reservoir tank. The UAM produces up to 10 liters per minute of up to 95% oxygen. The oxygen then travels through the flow meter into the back bar and reservoir bag.

GAS FLOW IN THE PATIENT SYSTEM

The flow meter controls the amount of oxygen entering the back bar and available for the patient. The spinning bobbin shows the flow rate in liters per minute.

The back bar provides the following functions:

1. The reservoir bag stores up to 2 liters of oxygen.
2. The pressure relief valve limits pressure in the back bar to 5 centimetres of water, protecting the patient from harm.
3. The air inlet valve allows room air to be drawn in whenever the supply of oxygen is less than the demand, as expressed by the patient's minute volume.

The low-resistance draw-over vaporizer adds volatile agent to the carrier gas passing through it. Volatile agent is added by using the selector wheel to achieve the desired percentage. The resulting anaesthetic supply gas is comprised of oxygen, volatile agent and a variable amount of room air. The anaesthetic supply gas flows from the back bar at low pressure to the breathing block and the bellows.

An oxygen monitor measures the oxygen concentration of the supply gas passing to the patient. The reading is displayed on the control screen.

As the supply gas flows from the oxygen sensor it passes the first one-way valve, then the bellows chamber, followed by the second one-way valve.

During spontaneous breathing both non-return valves open for inspiration and close during expiration. The bellows will not fill with exhaled gas.

The bellows is used for intermittent positive pressure ventilation, also known as IPPV or controlled ventilation. When the bellows is raised, supply gas flows through the first non-return valve into the bellows. The second non-return valve remains closed.

As the bellows are pushed down, the supply gas flows out through the second non-return valve into the inspiratory limb. The first non-return valve is closed.

Another pressure relief valve ensures that excessive manual force on the bellows cannot cause harm to the patient.

The UAM uses a conventional Y-piece patient connection. Exhaled gas returns to the UAM via the expiratory limb of the Y-piece and then passes by the Fenton balloon and a third one-way valve to the scavenger port. The balloon acts as an inflating valve which enables controlled ventilation of the lungs, the third valve prevents back-flow during spontaneous breathing.

SPECIFICATIONS

PHYSICAL AND ELECTRICAL

- Workstation: 146cm x 53cm x 69cm, 130kg, aluminium frame, vertical dove mount side rails on both sides, nylon internal tubing
- Top Shelf: 46cm x 31cm, 88cm above floor, 35kg weight limit
- Drawer: 10cm x 29cm x 32cm, 35kg weight limit, stainless steel, removable for cleaning.
- Top shelf: Epoxy powder coated pressed aluminium; side slots for nylon strap to secure equipment; optional monitor brackets are available to bolt items to the top of the top shelf, 54 cm (W) x 37 cm (D), Load maximum: 35 kg with even load distribution.
- Work surface: Stainless steel, removable for cleaning, 46 cm (W) x 31 cm (D) x 88 cm (H) above ground
- Mains power supply: 220V, 50-60Hz
- Power: 500 Watts/220v = 2.27 amps
- Control screen power: CE-marked, fused and medical grade. Mains powered with battery backup for one day of use.
- Automatic over/under voltage mains power isolator protects oxygen concentrator, and sockets for monitor and accessories
- 3.5" TFT touch screen for oxygen display and alarm setting
- Membrane switch for oxygen and alarm settings
- Casters: 150 mm diameter single wheel antistatic casters, front lockable casters

OXYGEN SUPPLY AND MONITORING

- Oxygen concentrator flow rate: 0.1 to 10.0 liters per minute
- Maximum oxygen concentration: 95%

NOTE: The oxygen concentration may vary according to ambient humidity and maintenance of the air filter

- Alternative sources: cylinder yoke, pipeline, or other external flow regulated portable oxygen source
- Automatic room air intake when patient minute volume exceeds supply gas flow
- Accuracy of glass rotameters: +/- 2.5% when using 100% O₂
- Integrated inspiratory oxygen monitoring uses MOX-3 oxygen sensor
- Calibration for room air (21%) and 100% oxygen
- Pressure sensitive apnea or high flow alarm

- Adjustable minimum and maximum oxygen alarm settings
- Up to ten hour battery backup, trickle recharge from mains power
- Membrane keypad and touch screen

OPTIONAL USE OF NITROUS OXIDE


- Sources: pipeline or cylinder
- Hypoxic cut-off: A solenoid automatically shuts off nitrous oxide delivery if supply gas O₂ level falls below a minimum of 25%
- Flow of N₂O stops if electrical power fails
- Sight and touch differentiated flow control knob per ISO standards
- Accuracy of glass rotameter for N₂O: +/- 2.5%

VENTILATION


- Bellows for adult and paediatric use made from durable, long-lasting silicone rubber
- 1600cc capacity
- Inspiratory pressure relief of 55 cmH₂O for pressure created during mechanical ventilation
- Aluminium bellows block
- Silicone balloon inflating valve at the expiratory port
- 15mm female/22mm male ISO standard taper connection for breathing circuit 30mm male ISO standard taper connection for AGSS scavenging equipment

VAPORIZERS


- Stainless steel and plated brass construction
- Draw over flow type
- Pour filler type.
- Separate units for isoflurane and halothane per ISO standards
- Agent delivery range: 0.5% to 4%
- 120ml capacity
- Performance:

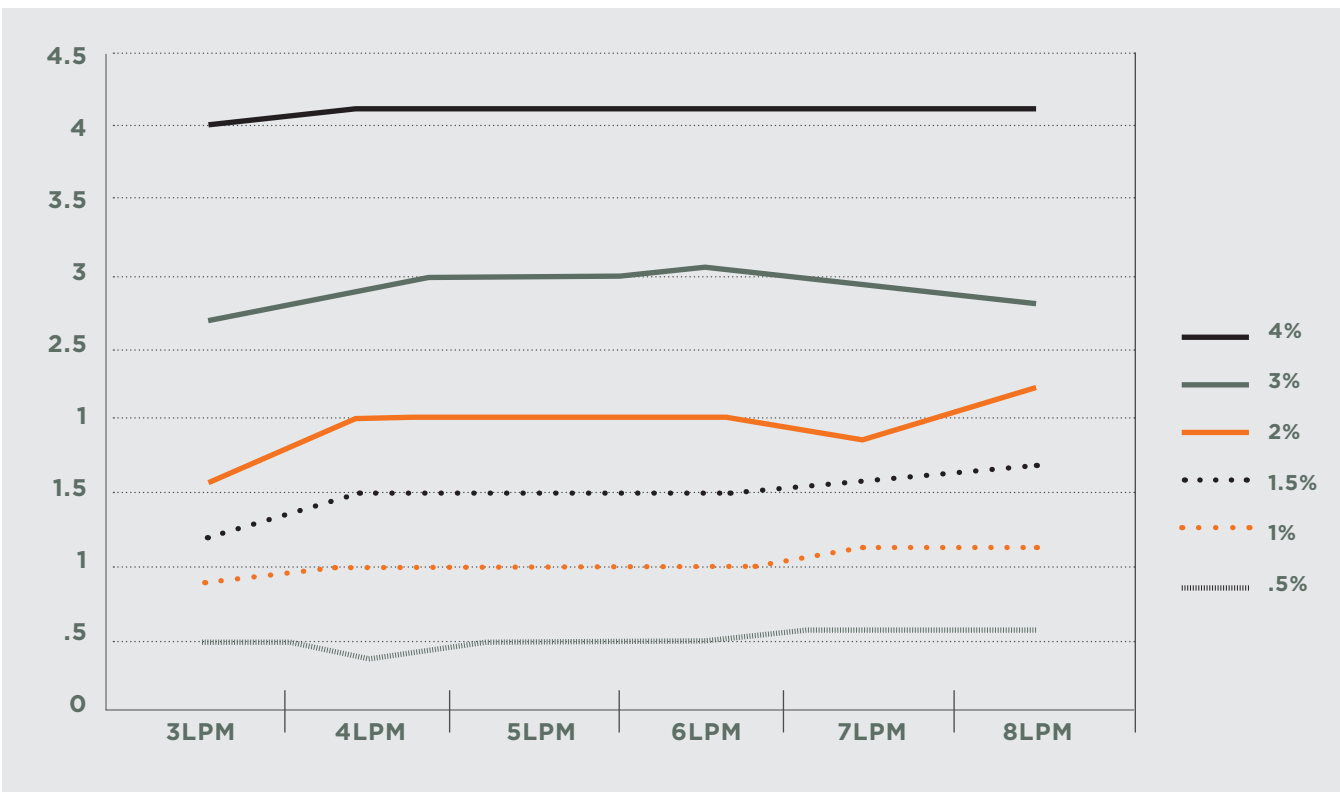


TYPICAL HALOTHANE VAPORIZER PERFORMANCE WITH FLOW

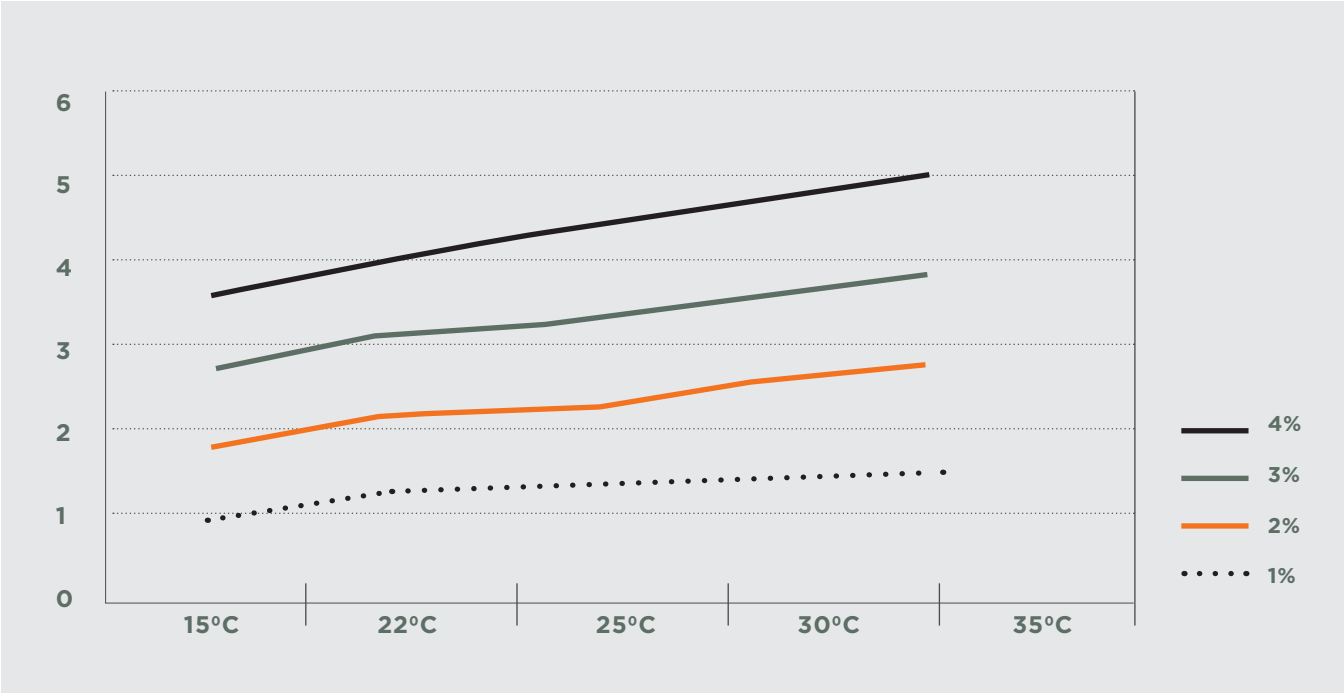

TYPICAL HALOTHANE VAPORIZER PERFORMANCE WITH TEMPERATURE

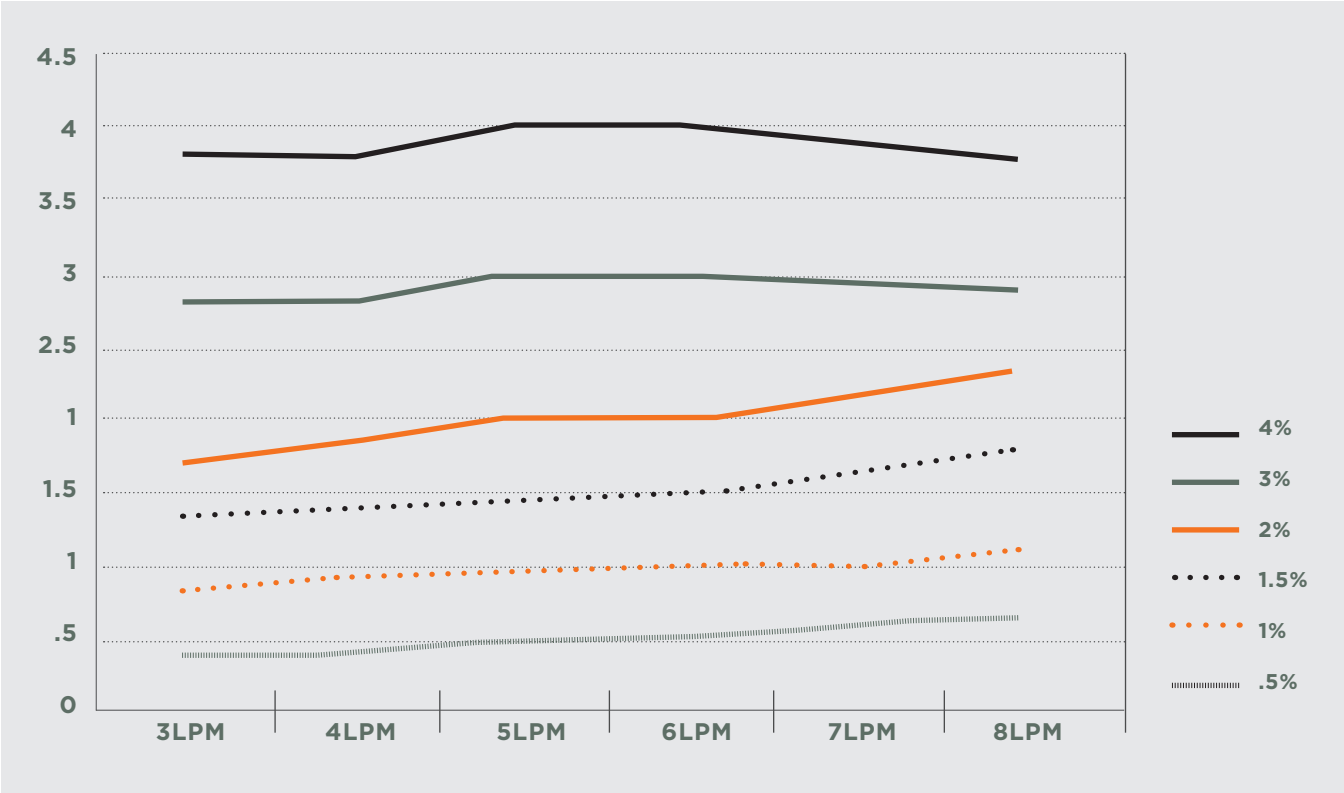

TYPICAL ISOFLURANE VAPORIZER PERFORMANCE WITH FLOW

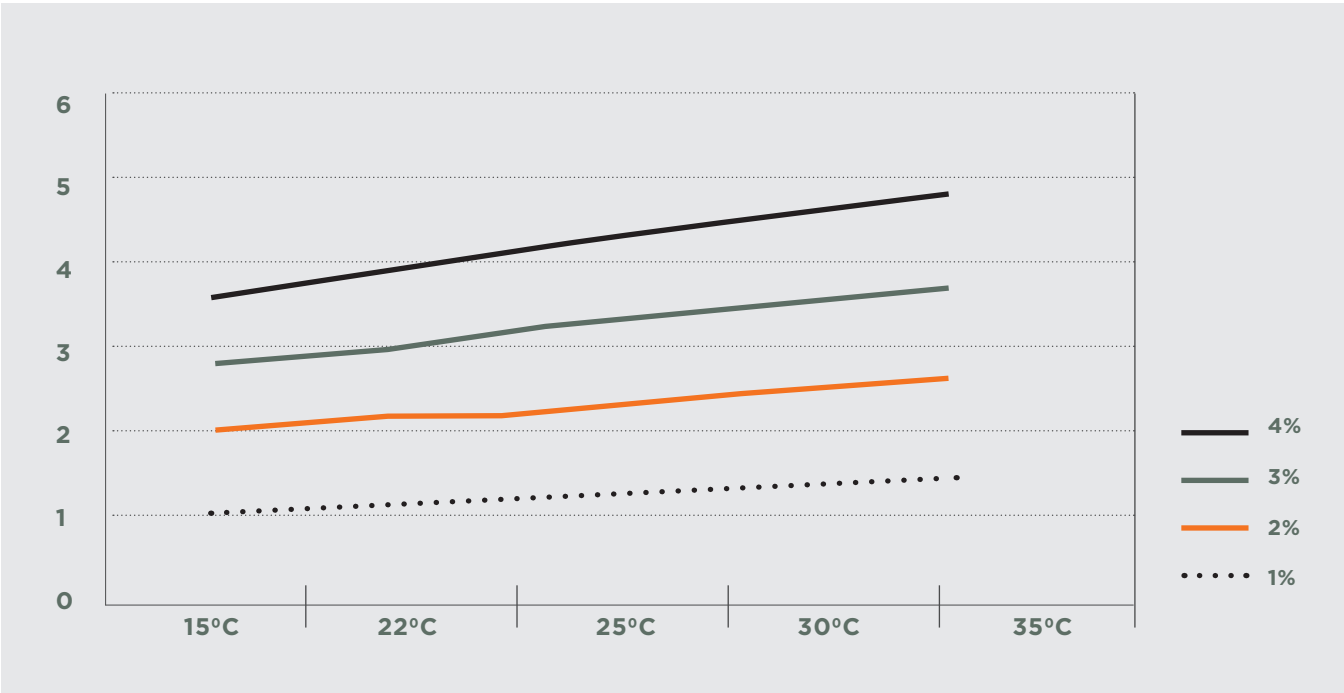

TYPICAL ISOFLURANE VAPORIZER PERFORMANCE WITH TEMPERATURE

REAR PANEL CONNECTORS


- Mains isolator switch, 2 UK style power sockets
- Oxygen and nitrous pipeline NIST connectors and cylinder yokes
- Multi-diameter external oxygen hose connector

COMPONENTS OF THE UAM



SEE INSET PAGE 7

INSET

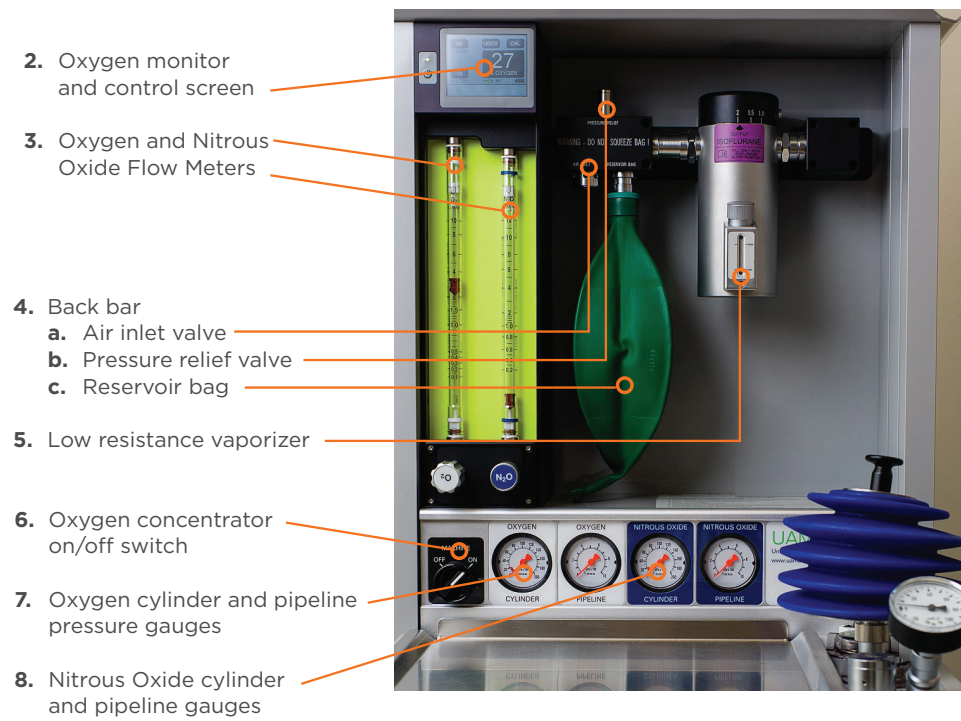

Oxygen sensor

Bellows

Pressure relief gauge

Breathing tubing connectors

Fenton Balloon
expiratory valve



Pressure relief valve

Breathing tubing connectors

Fenton Balloon
expiratory valve

Water trap





Fenton Balloon
expiratory valve

One Way Valve and
scavenging 30mm
taper connection

Percent Oxygen (shows the
oxygen percentage of the
anaesthetic gas measured
before inspiration)

THE CONTROL SCREEN

Calibration settings (to
calibrate oxygen sensor
and apnea alarm)

High oxygen percentage
alarm (user adjusted)

Mode (opens screen for
system settings)


Low oxygen percentage
alarm (user adjusted)

Battery charge status

Message Area

Apnoea alarm on/off

Mains power indicator
(machine is on mains power supply)

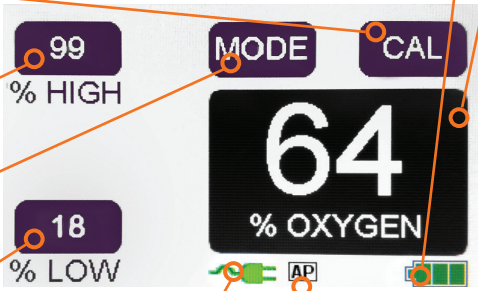

Pin-indexed connectors for E-size oxygen and nitrous oxide cylinders

Auxiliary oxygen inlet

Pipeline connections for oxygen and nitrous oxide

Mains isolator (Power On/Off)

Fuses for oxygen monitor, oxygen concentrator

British style fused 3-pin plugs



UNPACKING THE UAM AND PREPARING FOR USE

1. Remove front panel from the crate and remove the internal wooden crossbeam. If applicable remove the vital signs monitor box and set it aside.
2. Remove accessories from underneath and along the sides of the UAM and set aside.
3. Unlock the front wheels and carefully roll the machine out of the packing case, supporting the weight as you lower the machine to the floor. Remove all packing material.
4. Unpack work tray from its cardboard packaging. Remove film from work tray and place on top of the drawer.
5. Connect pipeline hoses and cylinders.
6. Screw in oxygen sensor and connect telephone-style cable. Some oxygen sensors come with an extra piece: this may be discarded. (See photo on next page.)

USING OXYGEN

The UAM accepts an E-size pin-index oxygen cylinder. The cylinder connects via a yoke and Bodok seal and the pressure is displayed on the pressure gauge on the front of the machine.

Bull-nose cylinders may be attached to the pipeline inlet using an approved 4 Bar medical regulator.

USING A PIPELINE CONNECTOR

The UAM accepts pipeline connections for oxygen. Each machine is supplied with an oxygen hose connecting to the UAM's NIST inlet. The other end is typically shipped with a BS 5682 probe, which is the British standard. Probes for other connector types can be arranged. Pipeline pressure is displayed on the gauge on the front of the machine.

USING THE AUXILIARY OXYGEN INLET

An auxiliary oxygen supply connection is mounted on the back of the machine. It accommodates differing tube sizes from ward-style flow-meters.

CAUTION: Oxygen from this source is not controlled by the glass flow meter on the front of the UAM. Rather, it must be controlled with an external flow meter attached to the cylinder providing the auxiliary oxygen. The user can confirm that oxygen is being given by observing the oxygen monitor.

NOTE: Use only 1 liter/min flow to conserve the oxygen in the external cylinder. This will provide 30-35% oxygen to the patient.

USING NITROUS OXIDE

The UAM supports cylinder and pipeline sources of nitrous oxide. The flow of N_2O is controlled by the right hand side flow meter. The N_2O flow knob has a different size and shape to distinguish it from the oxygen flow meter control knob.

The UAM automatically cuts off nitrous oxide if the percent of oxygen in the anaesthetic mix falls to 25% and below. The oxygen monitor must be powered on and working correctly for nitrous oxide delivery to occur. An exhausted oxygen sensor or battery will shut off N_2O supply.

The cut off mechanism is electrical, unlike older style mechanical interlocks. The N_2O valve is closed in its resting position and must receive a current from the oxygen monitor to open. Any interruption of that current, including a decrease in oxygen below 25%, will shut off the flow of nitrous oxide.

NOTE: Do not open the package of the spare oxygen sensor until it is needed. The oxygen sensor will start to deplete as soon as it comes into contact with air.

6. Screw in oxygen sensor and connect telephone-style cable. Some oxygen sensors come with an extra piece: this may be discarded.




7. Attach bellows assembly to UAM by placing it on top of the connector and firmly pushing down. Then tighten knurled nut fully.



8. Attach airway pressure gauge by pushing down collar and inserting gauge into the port. Press firmly until collar springs up.



9. Attach green reservoir bag to back bar as shown in photo.

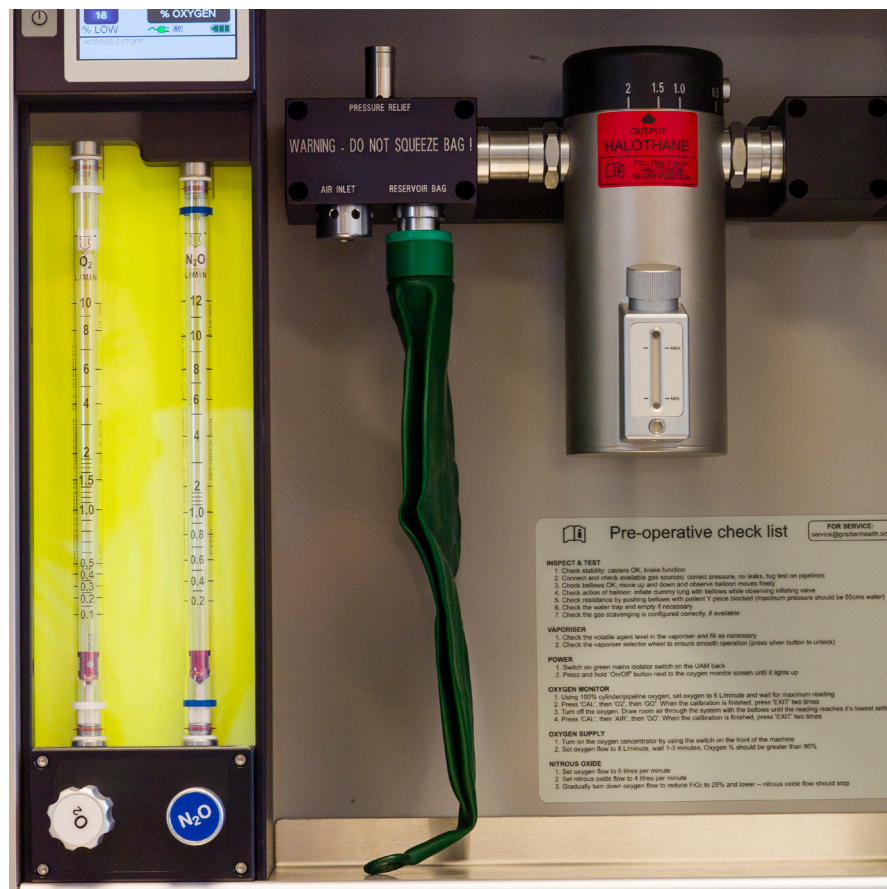


10. Attach patient circuit to inspiratory and expiratory ports.



11. Attach green plastic 30M-30F connector to the scavenging hose by inserting the male side into the hose. The female side attaches to the one-way expiratory port underneath the Fenton balloon.
12. Connect 30mm hose to scavenging port, or attach standard scavenging device. The use of a scavenging pressure relief valve is highly recommended to prevent pressure build-up due to unforeseen obstructions.
13. Operate the bellows, check that it fills freely with air from the air inlet and observe the Fenton balloon moves freely in it's housing when the bellows is moved. Occlude the Y-piece and check airway pressure reaches 45-55 cmH₂O.
14. Connect and turn on mains power by pressing the green switch on the back of the UAM. Wait about 10 seconds for the internal supply to come on, and then check the screen to see if the oxygen monitor shows green. This indicates the mains supply is working correctly.

NOTE: GRADIAN HEALTH SYSTEMS STRONGLY RECOMMENDS CONNECTING THE UAM TO A 2000-WATT VOLTAGE STABILIZER.

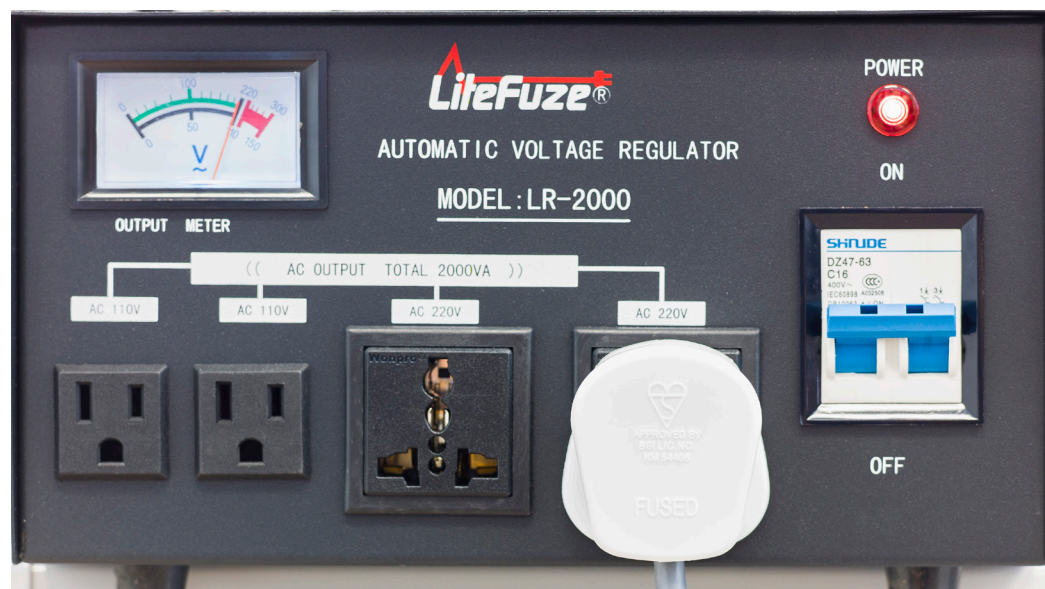

15. Calibrate oxygen monitor using 100% O₂ from cylinder or pipeline if available. NOTE: If cylinder or pipeline do not deliver 100% O₂, only calibrate on room air as described below in 15c.
- a. Using 100% cylinder/pipeline oxygen, set oxygen flowmeter to 6 L/minute and wait for maximum reading on monitor.
 - b. Press CAL, then O₂, then GO. When the calibration is finished, press EXIT twice.
 - c. Turn off the oxygen. Draw room air through the system with the bellows until the reading reaches its lowest value.
 - d. Press CAL, then AIR, then GO. When the calibration is finished, press EXIT twice.
16. Calibrate Apnea alarm pressure transducer as follows:
- a. Make sure flow meters are set to zero flow and reservoir bag is empty
 - b. Press CAL, then PRESSURE: CAL, then GO. When the calibration is finished, press EXIT twice.
17. Turn on oxygen concentrator using the black switch on the front of the machine. Set the oxygen output to 6 liters per minute and ensure that oxygen reaches at least 90 to 95% output after a few minutes.
18. If using nitrous oxide, check the function of the Nitrous Oxide cut off at 25% oxygen concentration as follows:
- a. Set oxygen flow to 6 liters per minute.
 - b. Set Nitrous Oxide flow to 4 liters per minute.
 - c. Gradually turn down the oxygen flow to reduce FiO₂ to 25% and lower—nitrous oxide flow should cut off.
19. Fill the vaporizer to maximum with correct anaesthetic agent. See the directions in the User Manual. **WARNING:** Do not overfill the vaporizer. Cover the oxygen sensor with a protective cloth. Do not spill agent on to the UAM, especially on to the oxygen sensor. Refer to the fill level window.
20. Place all manuals and training materials in drawer and secure a location for extras, such as the additional oxygen sensor, spare balloons, and so forth.

NOTE: High levels of relative humidity in the environment will affect the performance of the UAM. It is recommended that the UAM be used in areas where the humidity can be controlled. After operating the UAM in humid environments water drainage should be observed underneath the machine once the concentrator is turned off and has depressurized.

IF YOU HAVE PURCHASED THE PATIENT VITAL SIGNS MONITOR:

- 1. Unpack monitor and insert the battery. The battery is in a white box inside the foam protectors holding the monitor securely in the cardboard box.
- 2. Mounting the monitor: Please refer to the mounting instructions for the particular model of monitor obtained.
- 3. Connect patient cables to monitor – NIBP, SPO2, TEMP and ECG
- 4. Connect mains lead from monitor to one of the three prong sockets on the back of the UAM.
- 5. Allow batteries in patient monitor and control screen to charge for 24 hours.

RECOMMENDED UAM MAINTENANCE SCHEDULE

FREQUENCY	ACTION
Daily or prior to use (user)	Perform operational check (see UAM User Manual and Checklist on the UAM)
Every 6 months	Perform full function test Check air filter- clean or replace if needed
Every 5 years	Contact Grian Health Systems for a complete maintenance check service@gradianhealth.org +1.212.537.0340 +1.929.280.0210

FULL FUNCTION TEST (SEE CHECKLIST IN APPENDIX IV)

SYSTEM CHECKS

1. Check that no unauthorised modifications have been made to the UAM. Check that the Ayre's T-piece is present and spares as originally supplied are present. Check oxygen tubing for venturi (Hudson type) mask is available.
2. UAM frame – check for any loose screws or panels and tighten as required. Check wheels.
3. Check mains cable and gas supply tube condition and repair as required.
4. Connect mains supply, turn on mains isolator (note that there will be a delay of around 10 seconds while the automatic voltage switcher [AVS] starts) and check that the control screen and patient monitor (if fitted) yellow charging LED's are visible.
5. Turn off the machine mains isolator and ensure that the mains power Icon goes red – this indicates mains failure.
6. Check that the battery indicator icon shows three green segments and that it discharges with mains off and recharges with mains on. NOTE: Battery must be replaced every three years or sooner if battery life is low.
7. Check that the cylinder yoke seals are fitted and replace if damaged or missing.
8. Connect the cylinders to the yokes and ensure that the cylinder pressure is indicated on the correct pressure gauge (nitrous oxide and oxygen) when turned on individually.
9. Connect the machine pipeline hoses to the supply outlets (wall outlet or cylinder regulator) and ensure that pipeline pressure is indicated on the correct pressure gauge (nitrous oxide and oxygen) when connected and turned on individually.
10. Turn on the oxygen and the nitrous oxide flowmeter control knobs using each gas supply (concentrator, cylinder and pipeline) and ensure that the maximum flow of 10LPM is achieved for each gas.
11. Close the flowmeter control knobs and ensure that the flow stops.
12. Check the oxygen concentrator air filter and clean or replace if required.
13. Check the cooler fins and clean if required.
14. Test the oxygen concentrator loss of power audio alarm and replace with 9 V battery if necessary.

OXYGEN CALIBRATION, SAFETY AND QUALITY CHECKS

15. Turn on the control screen and ensure that the vaporizer is off. Turn on the oxygen cylinder or pipeline to provide 100% oxygen, turn on a 6 liter flow and wait for the oxygen sensor to stabilize, then press CAL – Oxygen – GO then wait, then back to the main display. If the oxygen sensor fails to calibrate or an oxygen sensor exhausted alarm is present replace the oxygen sensor.

NOTE: For the most accurate performance over the entire range of possible oxygen concentrations, the oxygen sensor must be calibrated with both 100 % oxygen and 21% oxygen (room air). If 100% oxygen is not available then the cell should at least be calibrated with 21% oxygen. If the sensor has not been calibrated at 100% oxygen then the monitor will display 104% at higher oxygen concentrations. Calibrating the cell at 100% oxygen can clear this condition.

NOTE: The oxygen sensor will need to be replaced approximately every 12 months. An alert on the screen will indicate when the oxygen sensor is exhausted.

16. With the oxygen turned off cycle the bellows to draw in air, wait for the oxygen sensor to stabilize, then press CAL – AIR – GO then wait, then press EXIT to go to the main display. If the sensor fails to calibrate and a oxygen sensor exhausted alarm is present replace the oxygen sensor.
17. Check that the upper and lower alarms can be adjusted and that they function, this is achieved by turning up the oxygen supply to exceed the high alarm or increasing the entrained air to fall below the lower alarm.
18. Increase the nitrous oxide flow, if available, and ensure that it cuts off when the oxygen concentration is 25% or lower.
19. Disconnect the oxygen sensor (mounted behind the bellows) – ensure that the disconnect alarm sounds and that nitrous oxide (if fitted) turns off.
20. Set a 6 LPM flow from the oxygen concentrator and ensure that the output stabilizes at 90% – 96% oxygen concentration, increase the flow rate to 10 LPM and ensure that the output remains at 90% or better. If the output concentration is low check the filter condition and clean or replace it if required.

NOTE: The above verification can only be performed if the oxygen sensor has been calibrated with 100% oxygen.

BREATHING SYSTEM CHECKS

21. Turn on the oxygen flow with patient Y-Piece blocked and ensure that the reservoir bag fills and excess gas escapes through the pressure relief valve and that maximum circuit pressure is around 5 cmH₂O. Confirm that the apnea alarm sounds with a full bellows after 30 seconds and that it repeats every 30 seconds.

22. Turn off the flowmeter and cycle the bellows and ensure that air is drawn in through the under pressure valve - this is indicated by a slight fluttering of the reservoir bag when it is empty.
23. Check condition of reservoir bag.
24. Check that the bellows assembly is securely connected to the machine and is undamaged, check for splits in the bellows or displacement of the end plates, loosen and re clamp centrally if required.
25. Cycle the bellows assembly with the inspiratory port blocked and ensure that 55 cmH₂O can be achieved when the bellows is pushed down for inspiration.
26. Check the patient system water trap is correctly screwed in place without leaks.
27. Check the correct function of the Fenton balloon. Attach the Y-piece patient connection and occlude the distal end. Operate the bellows and observe the free movement of the balloon in its tube. Apply sustained pressure on the bellows and check the pressure rises to at least 45 cm H₂O. If patient system pressure is not achieved, unscrew the balloon cover, remove balloon and perform checks (see page 20)
28. Cycle the bellows and ensure that the pressure relieves at 55 cmH₂O.

VAPORIZER CHECKS

29. Check that the vaporizer is full, fully turned off and that it indicates the level of liquid.
30. Turn on the concentrator and set an oxygen flow of 6 liters per minute and wait until the output has stabilised. Set the vaporizer output to 3%. Verify that after 2 minutes that the oxygen concentration has dropped by 3%.

NOTE: This is a very basic test and can only indicate approximate output. The vaporizer output is affected by temperature, flow and time. Ideally, the above test should be conducted at an ambient temperature of 22 degrees Celsius.

RESIDUAL FLOW CHECK

31. Turn on the concentrator and allow it to run for at least 5 minutes with the oxygen rotameter fully closed. Turn the concentrator off and set the oxygen rotameter flow to 200 ml/min. You should observe this level of flow for at least 5 minutes. If this test fails this could indicate a leak anywhere between the oxygen concentrator reservoir tank and the rotameter.
32. Ensure that the service record (see appendix IV) is completed and filed. Attach a label (see appendix V) with your initials and the date to let the users know that the UAM has received periodic maintenance and testing.
33. Perform ventilator checks if available as indicated in the ventilator maintenance manual and separate checklist.
34. Perform patient monitor checks if available per separate monitor checklist.

MAINTENANCE PROCEDURES

CLEANING OR REPLACING THE OXYGEN CONCENTRATOR FILTER

Materials and tools required:

- 3mm hexagonal wrench
- M5 flare nut wrench

The oxygen concentrator includes a replaceable filter that cleans the room air entering the oxygen concentrator. Keeping it clean will extend the life of the concentrator.

To remove the filter:

1. Turn off the UAM and disconnect it from the mains power supply.
2. At the lower back panel, loosen the six hex screws using a 3mm hexagonal wrench.
3. Carefully tilt the rear panel back from the top, taking care not to disconnect the wires attached to it.
4. Detach the two hoses leading to the air filter. Press the orange collar to release the hose. Do not pull on the hose without depressing the collar.
5. Remove the cover plate that holds the filter in place by removing the centre M5 nut.
6. Filters may be cleaned by shaking or vacuuming. If the filter is too dirty to clean, replace it with a new one.
7. Install the new filter and replace the cover plate.
8. Inspect the air cooler fins located in between the sieve tanks. Remove any dust with a damp cloth.
9. Replace the back panel, connect the mains power supply, and turn the UAM on.

NOTE: Replacement filters can be ordered from Gradian (Part number: 1300-069) but they may be less expensive and quicker to obtain from local automotive parts suppliers.

Filter Details:

Type: AG285

Applications: Used in automobiles made by Ford, Audi, Volkswagen, Citroen, Seat, Skoda and Saab.

Replacement options:

Factory Number

FRAM	CA4739
GMC	93152533
MEYLE	7126080390
MONARK	30780031
PUROLATOR	AF3561
TEHO	2247

Dimensions:

Outer Diameter: 241 mm

Inner Diameter: 185 mm

Height: 50 mm

REPLACING THE OXYGEN SENSOR

Materials and tools required:

→ New sensor

Replace the oxygen sensor every 12 months, or when an alarm message on the oxygen sensor indicates that the sensor has failed to calibrate.

Always keep a spare sensor and store it unopened in its sealed pouch. The sensor begins to decay when you open the packaging and expose it to air.

To replace the oxygen sensor:

1. Disconnect the oxygen sensor cable.
2. Turn the existing sensor anticlockwise to disconnect it.
3. Remove the new sensor from its container and remove the extension piece covering the bottom of the sensor.
4. Place the sensor in position, and turn it clockwise to screw it into place. Do not over tighten.

NOTE: Take care to engage the thread correctly, and ensure that the replacement cell is fitted with a suitable O-ring to prevent leaks from the breathing system. Replacement sensors include a new O-ring.

5. Connect the oxygen sensor cable.
6. Recalibrate the oxygen sensor.


REMOVING AND REPLACING THE DRAW-OVER VAPORIZER

CAUTION: Do not use solvents to clean the oxygen sensor.

Materials and tools required:

→ 4mm Allen key

1. Ensure that the UAM is turned off.
2. With a fine tipped pencil draw a line on the back bar along the right edge of the right hand block. This will help align the block when re-assembling.
3. Use a 4mm Allen key to turn the right hand block's 4 screws counter clockwise.

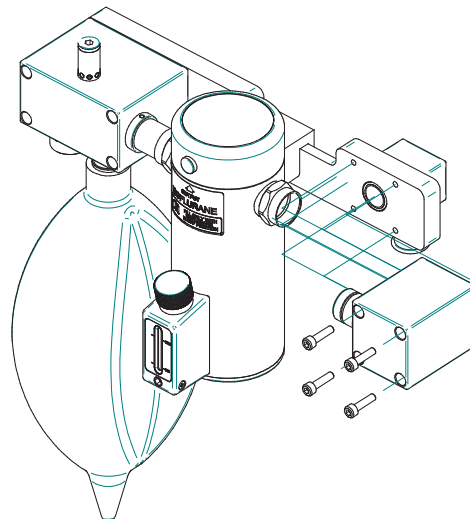

4. Slide the block off the vaporizer.
5. Slide the vaporizer to the right to remove it from the back bar.
6. Lift the vaporizer off the back bar.
7. To replace the draw over vaporizer: repeat the above steps in reverse order using the pencil mark as a guide to align the right hand block.

NOTE: Ensure that the O-ring in the back of the right hand block is in position during reassembly.

CLEANING THE HALOTHANE VAPORIZER

Materials and tools required:

- Metal bowl or container

Drain Halothane vaporizers every 3 months to prevent thymol build up. Excess thymol can damage the selector wheel mechanism over time.

1. Unscrew the filler cap on the front of the vaporizer.
2. Place a metal bowl or container underneath the vaporizer. Halothane is corrosive and can damage the UAM if spilled.
3. Using the filler cap, insert the key on the bottom of the cap into the receptacle just below the window on the front of the vaporizer.
4. Turn the filler cap anticlockwise to open the drain underneath the window.
5. When finished, turn the filler cap clockwise until the drain is fully closed. Do not over tighten.
6. Re-fill the vaporizer with halothane.
7. Dispose of the old anaesthetic agent using your hospital's procedures for disposing of medical waste.



REPLACING THE FENTON BALLOON

Materials and tools required:

→ New Balloon

1. Unscrew the tube assembly unit, which includes the top aluminium ring, the clear plastic housing and the valve assembly attached to the bottom. Turn using the top ring to make sure that it comes off as one piece. Make sure that the O-ring above does not fall down.
2. With fingers on the neck of the balloon, carefully peel off the balloon from its mounting.
3. With a small amount of moisture on the neck, slide the new balloon on to its mounting and move it up so that the top ring locates in the groove on the mounting.
4. Replace the tube assembly unit, taking care not to catch the balloon in the screw threads. Use a rotating action to make sure the balloon is fully visible inside the tube before locating the screw threads.




REPLACING THE CONTROL SCREEN CIRCUIT BOARD

Materials and tools required:

- 4 mm Allen key
- 2.5 mm Allen key
- Long nose pliers

1. Switch off mains power.
2. Remove the six monitor shelf retaining screws and washers using a 4 mm Allen key and lift off the monitor shelf.
3. Disconnect the oxygen sensor cable - note that the clip must be squeezed in to disconnect.
4. Pull off the silicon tube from the pressure transducer - do not put any stress on the transducer - support the body and use a slight twisting action.
5. Disconnect the 4 way power supply connector - note that the clip must be squeezed in to disconnect.
6. Disconnect the speaker connector - pull connector body out with long nose pliers - do not pull on cable.
7. Disconnect the solenoid connector - pull connector body out with long nose pliers - do not pull on cable.
8. Remove the two screen securing screws using a 2.5 mm Allen key.
9. Slide screen up to remove from flow meter.
10. Pull touch screen connector clamp to disconnect the ribbon cable.
11. Disconnect screen display cable by flipping up the clamp using a thumbnail.
12. Remove the 4 printed circuit board securing screws using a 2.5 mm Allen key.
13. Reverse sequence using new board.

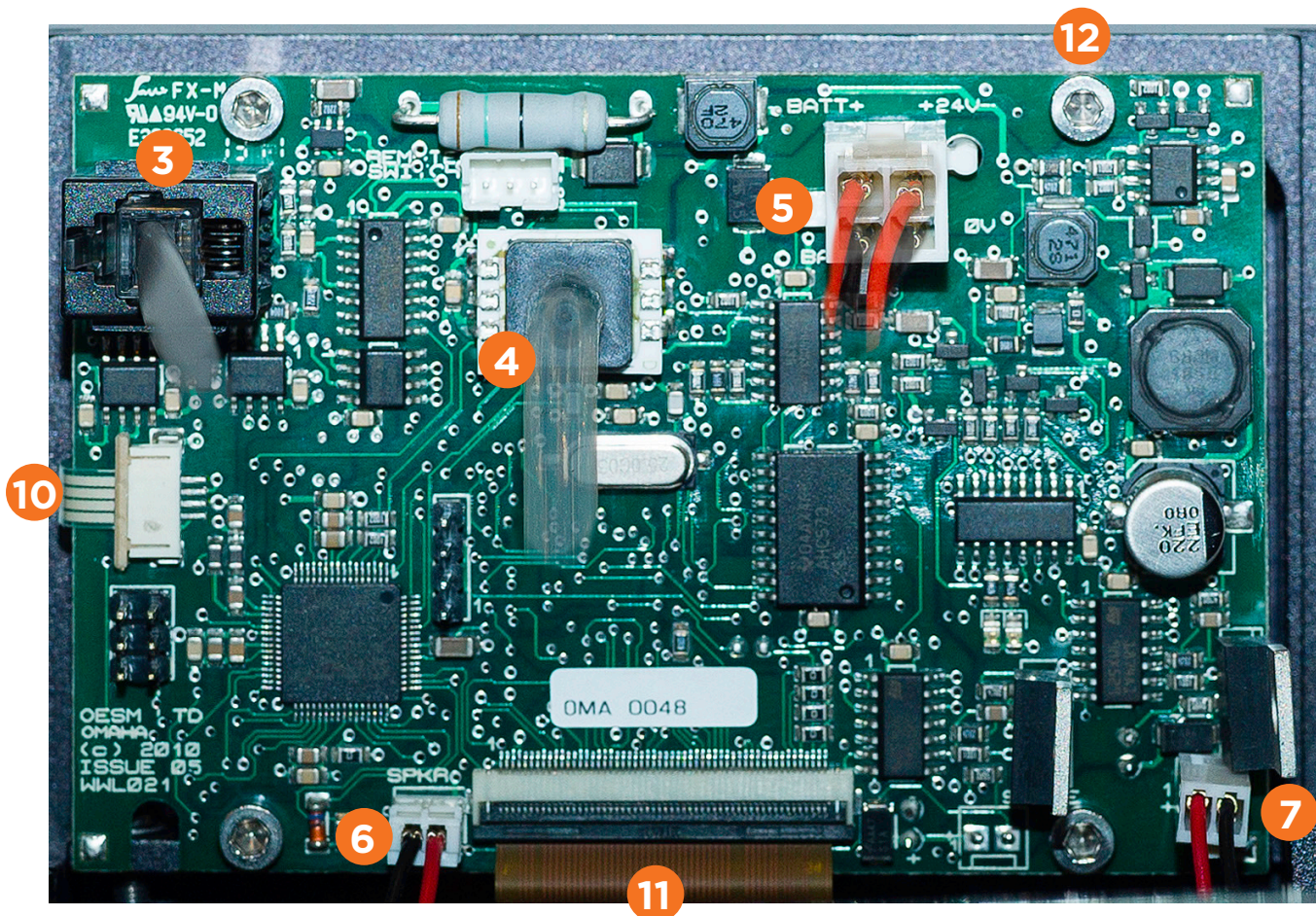


REPLACING THE CONCENTRATOR LOSS OF POWER ALARM BATTERY

Materials and tools required:

- 3 mm Allen key
- 9 volt battery

1. Turn off the UAM and disconnect it from the mains power supply.
2. Loosen the hex screws on the back lower panel.
3. Lift and remove the panel and place it so that it rests on the left side of the UAM making sure not to pull on the cables.
4. At the back left side of the concentrator, disconnect the battery
5. Install a new 9 volt battery. Reuse the Velcro strip from the old battery.
6. Replace the panel.

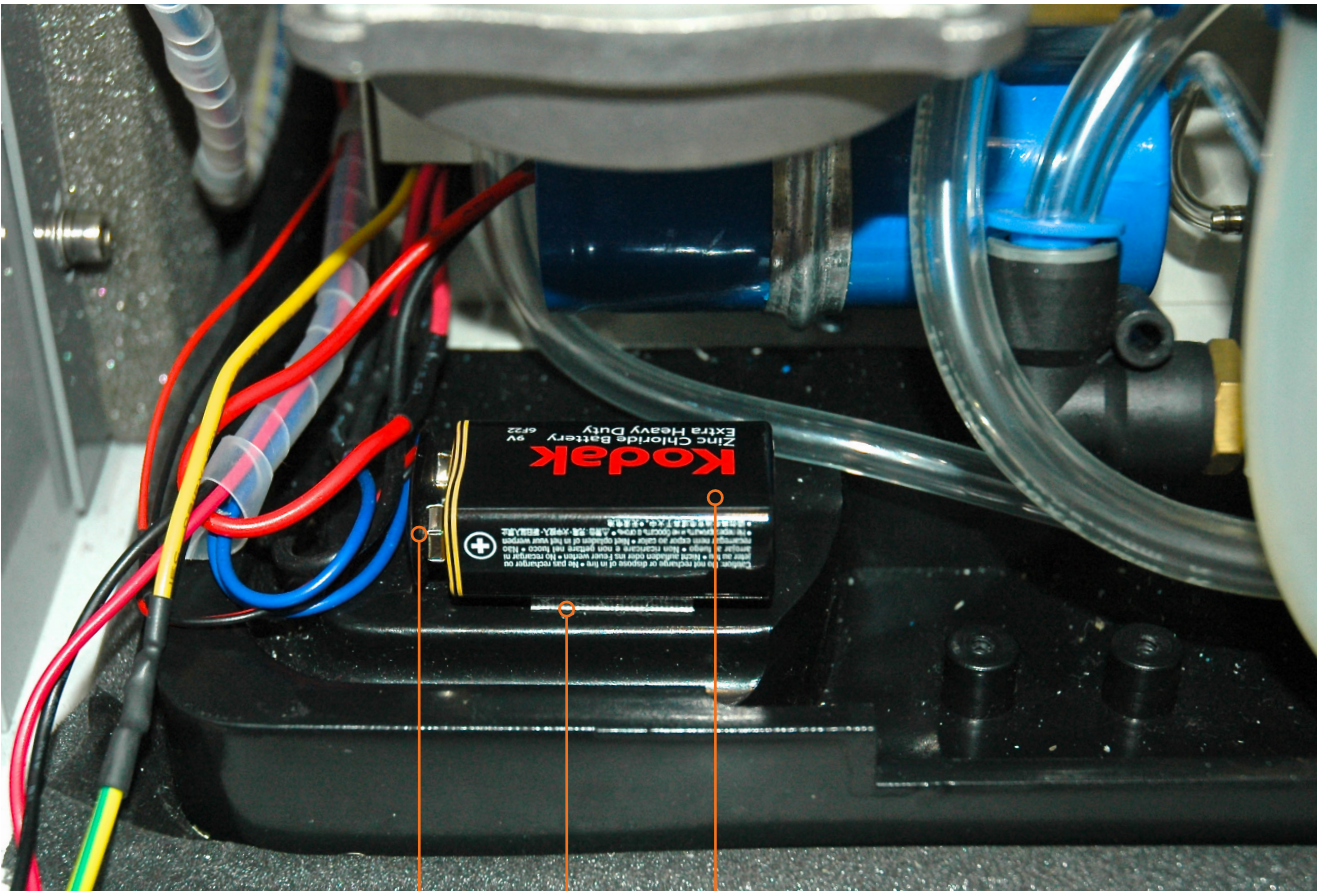

BATTERY CONNECTOR

VELCRO

9V BATTERY

REMOVING THE CONCENTRATOR

Materials and tools required:


- 6 mm Allen key
- 2.5 mm Allen key
- M5 flare nut wrench
- Long flat screw driver or other long and flat metal tool

1. Engage the brakes on both of the UAM's front castors.
2. Remove the six screws from the lower rear panel.
3. Tilt the panel back and remove both hoses from the air filter by carefully pushing the orange plastic rings towards the filter and pulling the hoses away.
4. Rest the panel against the left side of the UAM being careful not to stress the cables connected to the panel.
5. With the nut wrench unscrew the nut on the ground post terminal on the bottom right of the panel to remove the ground wire connecting the panel to the concentrator.
6. Disconnect the concentrator power switch cable (4 wires: 2 blue, 1 yellow, 1 red).




7. Disconnect the white nylon tube going from the concentrator to the bottom of the regulator by carefully pushing up the orange plastic ring from the metal connector and pulling the tube down.
8. Remove the (4) 6mm hex screws surrounding the rectangular opening underneath the UAM.
9. Push the top of the concentrator module tilting it towards the front of the UAM and with a very large flat screwdriver or another hard flat tool separate the bottom of the concentrator from the bottom surface of the UAM's interior.
10. Once the concentrator is separated from the base, lift it and pull it towards you and letting it rest on the back edge of the bottom surface of the UAM being careful not to pull on the concentrator monitor and LED indicator panel cables.
11. Disconnect the concentrator monitor cable and LED indicator panel cable.
12. Lift and lower the concentrator module to the floor right behind the UAM.

7

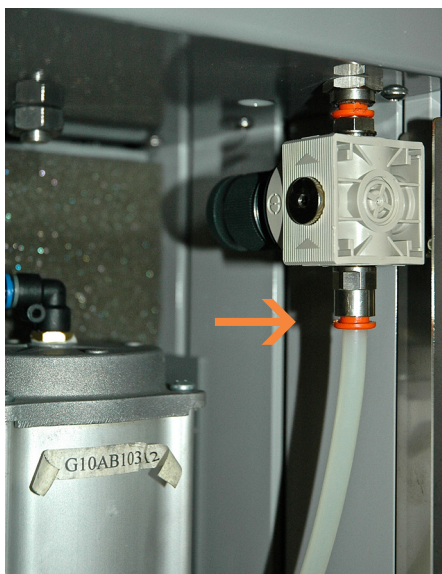


8

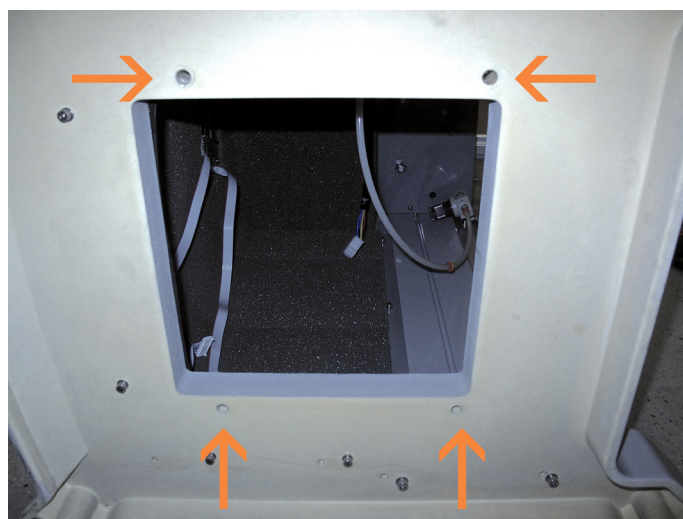

9



11

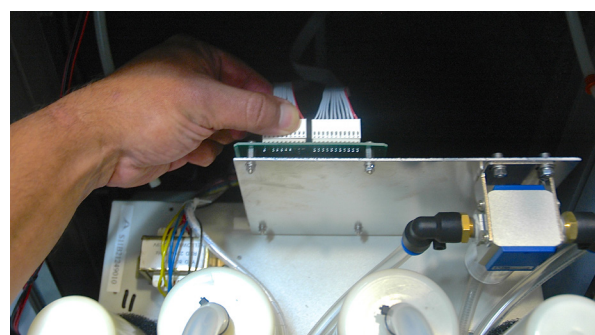

INSTALLING THE CONCENTRATOR

Materials and tools required:

- 6 mm Allen key
- 3 mm Allen key
- M5 flare nut wrench

1. Lift up and balance the concentrator on the rim of the chassis.
2. Reconnect the concentrator monitor cable and the LED indicator panel cable.
3. The drain hose from the water trap must be pushed through the hole in the bottom of the concentrator.
4. Push the concentrator into the housing, making sure that the drain hose does not get trapped between the bottom of the concentrator and the rim of the housing.
5. Make sure that the oxygen outlet tube from the oxygen reservoir is not kinked. Also make sure that the water trap is not damaged or dislodged by the top housing or by the foam on the right hand side of the housing. Realign it if necessary after it is inside.
6. When the concentrator is placed approximately in its final position, align the rear right hand holes for the 6 mm hex screw with a short screwdriver and loosely screw it in a few turns.
7. Align the rear left hand holes and loosely screw the second hex screw here. Then locate the remaining two front hex screws and tighten all screws.
8. Connect the oxygen supply hose to the bottom of the pressure regulator making sure that it is not kinked.
9. Connect the concentrator power switch cable (4 wires: 2 blue, 1 yellow, 1 red).
10. Make sure that the battery wires are connected firmly to the battery on the rear panel (black to black, red to red)
11. Insert the two air inlet hoses into the air filter and lift up the rear panel ensuring that the hoses bend downwards between the zeolite canisters and do not kink when closing the rear panel.
12. Attach the rear panel to the housing with the 6 screws.

NOTE: When finished assembling verify that the drain hose is protruding from the orifice as picture below.

13. Test the oxygen flow and concentration of the concentrator and test that the control screen works on battery and that the battery charges.



REPLACING THE CONTROL SCREEN BATTERY


Materials and tools required:

- 2.5 mm Allen key
- Long nose pliers

1. Remove the six screws from the lower rear panel.
2. Tilt the panel back and remove both tubes from the air filter.
3. Rest the panel against the left side of the UAM being careful not to stress the cables connected to the panel.
4. With the long nose pliers disconnect the positive and negative battery wires.
5. Remove the four screws on the back of the panel that attach the battery holder.
6. Remove battery and replace with a new one

NOTE: Replacement batteries can be ordered from Gradian (Part number: 1419-005) but they may be less expensive and quicker to obtain from a local battery supplier. Battery details: YUASA NP1.2-12 12V 1.2 AH sealed lead acid battery.

7. Assemble by following the above steps in reverse.

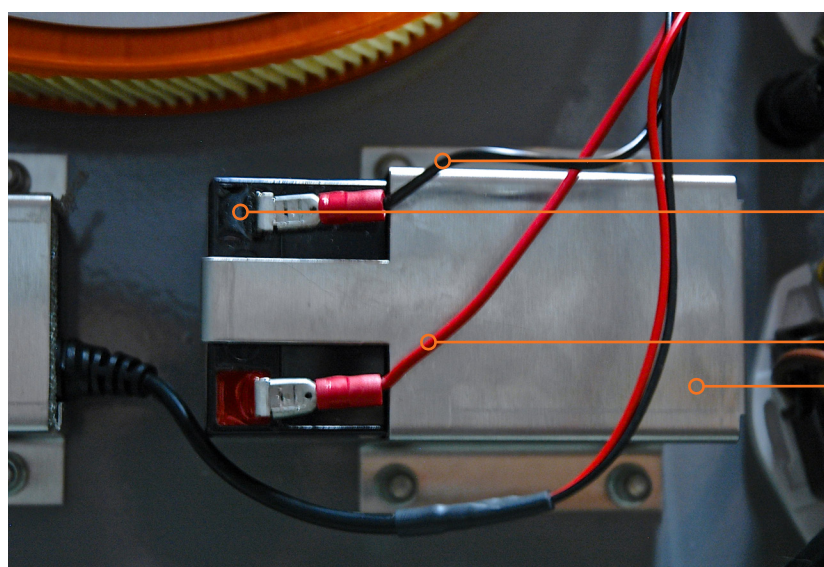


NEGATIVE BATTERY WIRE

CONTROL SCREEN BATTERY

POSITIVE BATTERY WIRE

BATTERY HOLDER

4

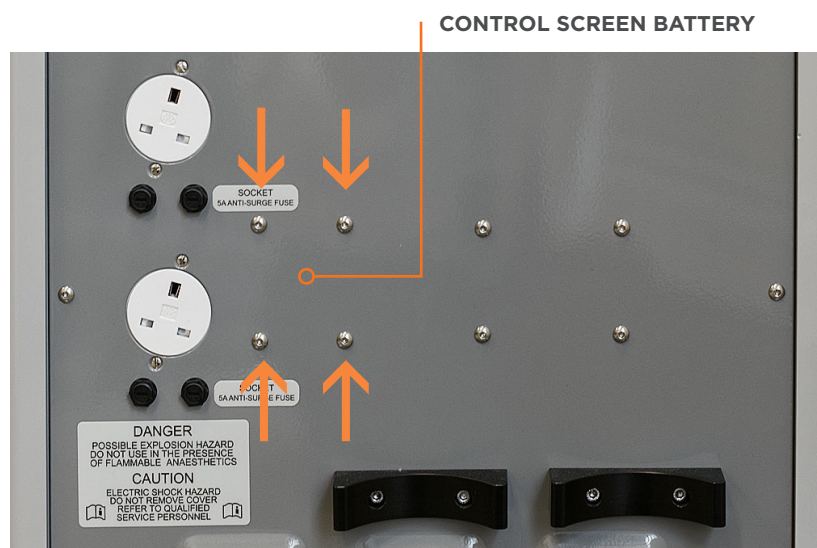

CONTROL SCREEN BATTERY

5



REMOVING, INSPECTING AND CLEANING THE PRESSURE RELIEF VALVE

Materials and tools required:

- 5mm hexagonal wrench
- A mole wrench (also known as vise-grips or locking pliers)
- A small cloth
- Medical grade cotton swabs
- Isopropyl alcohol (also known as rubbing alcohol or surgical spirit- typically 70% isopropyl and 30% distilled water)
- A small bowl

On occasion the 5 cmH₂O pressure relief valve located on the top of the back bar may not perform optimally due to accumulation of dust and other contaminants such as liquids. If you notice that the valve is not opening at 5 cmH₂O then you will need to remove it, inspect it, and clean it by following these steps:

REMOVING THE PRESSURE RELIEF VALVE:


1. Turn off the UAM.
2. Using a 5mm hexagonal wrench, loosen the plug on the top of the pressure relief valve and remove it.




NOTE: It may be possible that when removing the plug, the body of the pressure relief valve may loosen up from the back bar. This is not a problem and in fact eliminates the need for the next step.

3. With a mole wrench and a cloth, grab the body of the pressure relief valve and turn it counter clockwise to loosen it up and free it from the back bar.

NOTE: If you do not cover the pressure relief valve with a cloth before you grab it with the mole wrench, it may get damaged.

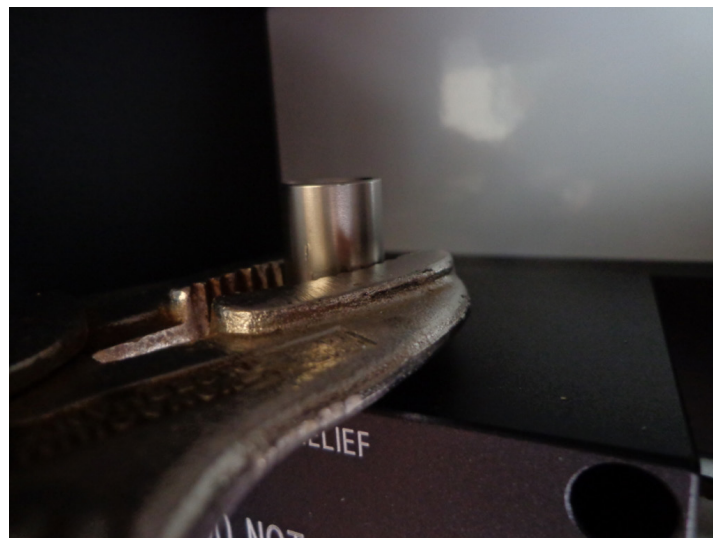

4. Once removed make sure that you see the following parts: the relief valve plug, the relief valve body, the stainless steel ball, and the O-ring attached to the valve's body.

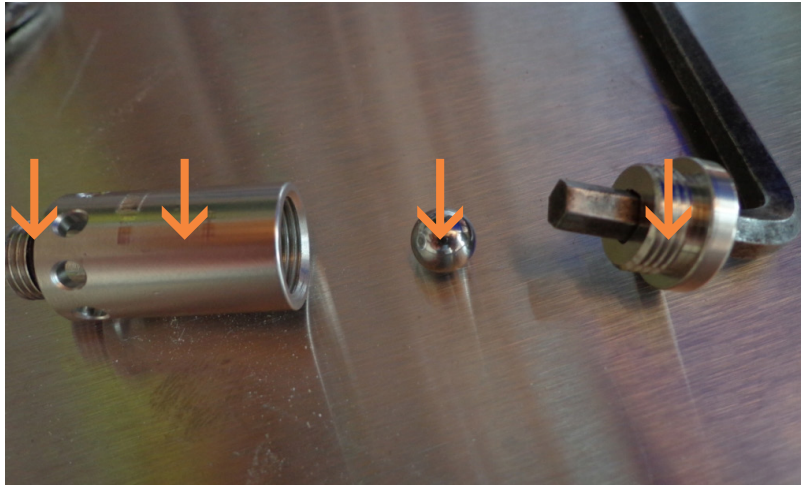

INSPECTING AND CLEANING THE PRESSURE RELIEF VALVE:

5. Inspect all of the parts to ensure that there are no stains, debris, moisture or rust. Inspect the condition of the O-ring and make sure it is not cracked or deformed.
6. To clean, place all the parts in a small bowl with isopropyl alcohol and let them soak for a few minutes.



7. Wipe each part with a medical grade cotton swab ensuring that there are no visible stains or debris.



8. Once cleaned, completely dry all of the parts.

REPLACING THE CLEANED PRESSURE RELIEF VALVE:

9. Place the stainless steel ball into the body of the valve and screw on by hand the plug onto the body.
10. With the 5mm hexagonal wrench screw in the assembled pressure relief valve into the back bar until tight, making sure not to exert too much force.
11. Test the operation of the pressure relief valve ensuring that it relieves pressure at 5 cmH₂O.
12. If the pressure relief valve does not function correctly, replace it with a new one.



PROBLEM-SOLVING

PROBLEM	POSSIBLE CAUSES/SOLUTIONS
Both the Oxygen concentrator and the control Screen do not turn on	<ul style="list-style-type: none"> → Check that both mains isolator switch and O₂ concentrator are switched on. → General supply voltage instability, power surges/low voltage/spikes. Use voltage stabilizer → Faulty socket outlet, general power cut → Malfunctioning mains isolator switch → Malfunctioning voltage switcher → Broken mains plug or cable → Malfunctioning voltage stabilizer, if fitted between UAM and wall socket.
UAM shuts down and re-starts.	<ul style="list-style-type: none"> → General supply voltage instability, power surges/low voltage/spikes. Use voltage stabilizer → Fault in local wiring or other device in or near the operating theatre. Check other electrical devices in the operating theatre such as air conditioner, sterilizer, theatre lights. Switch off all devices except UAM and switch on other devices at intervals, one by one.
Oxygen concentrator does not turn on or a noticeable hot smell	<ul style="list-style-type: none"> → Blown 5A fuse → Malfunctioning front switch → Malfunctioning mains isolator switch → Malfunctioning voltage switcher → Malfunctioning oxygen concentrator circuit board → Compressor failure
Control screen does not turn on	<ul style="list-style-type: none"> → Blown 24V Power Supply Unit fuse (500mA) → Defective 24V Power Supply Unit → Circuit board failure
Oxygen monitor displays 104% at higher oxygen concentrations	<ul style="list-style-type: none"> → Oxygen sensor needs to be calibrated with 100% oxygen
Oxygen sensor will not calibrate or monitor readings are erratic	<ul style="list-style-type: none"> → Oxygen sensor is exhausted and needs to be replaced → Sensor cable is disconnected or broken
Audible alarms not working	<ul style="list-style-type: none"> → Disconnected speaker connector → Malfunctioning speaker → Malfunctioning control screen circuit board
Control screen shuts off when mains power is lost	<ul style="list-style-type: none"> → Battery needs recharging or is disconnected → Battery has lost its capacity to charge and should be replaced
N ₂ O does not cut-off when oxygen percentage is less than 25%	<ul style="list-style-type: none"> → Malfunctioning cut-off solenoid valve → Malfunctioning control screen circuit board
Apnoea alarm does not activate	<ul style="list-style-type: none"> → Mute option is selected → Pressure transducer needs to be calibrated → Malfunctioning control screen circuit board
Oxygen concentrator output measures less than 90%	<ul style="list-style-type: none"> → Oxygen sensor needs to be calibrated to both 100% and 21% oxygen → Malfunctioning water trap in oxygen concentrator is allowing moisture to saturate the zeolite material in the Adsorption Towers → Air filter is dirty and obstructing air flow → Air compressor is malfunctioning → Switching valve is malfunctioning → Water in system due to clogged or defective water trap. Clean or replace. → Air filter is dirty and obstructing air flow. Clean or replace. → Air compressor is malfunctioning → Switching valve is malfunctioning
Oxygen concentrator output does not reach 10 liters/minute and/or flow fluctuates	<ul style="list-style-type: none"> → Concentrator output regulator needs to be adjusted
Concentrator alarm does not activate when power is lost	<ul style="list-style-type: none"> → 9-volt alarm battery is discharged and should be replaced.
Breathing system pressure too high	<ul style="list-style-type: none"> → Pressure relief valve not opening. Remove and clean. → Balloon valve malfunctioning. Remove and clean or replace.

MAIN PARTS

ITEM NO	DESCRIPTION	PART NUMBER
---------	-------------	-------------

GAS SUPPLY

1	Oxygen concentrator	8150-181
1-a	Main concentrator control board	N/A
1-b	Compressor pump	N/A
1-c	O ₂ concentrator switching valve	N/A
1-d	Zeolite Cylinder Set	N/A
1-e	Fan	N/A
1-f	Water Trap	N/A
2	O ₂ regulator assembly with yoke and contents gauge	8150-135
2-a	O ₂ cylinder pressure gauge	1515-004
2-b	O ₂ pipeline pressure gauge	1515-003
3-1	N ₂ O regulator assembly with yoke and contents gauge	8150-136
3-a	N ₂ O cylinder pressure gauge	1515-004
3-b	N ₂ O pipeline pressure gauge	1515-003
4	O ₂ 4 meter antistatic hose assembly with BS5682 connector	9950-017
5	N ₂ O 4 meter antistatic hose assembly with BS5682 connector	9950-027
6	Air filter- Type: AG2865	1300-069

ANAESTHESIA DELIVERY/BREATHING

8	Flowmeter clear screen	3300-225
9	Flowmeter knob O ₂ and label	2300-012 + 5750-002
10	Flowmeter knob N ₂ O and label	2300-013 + 5750-004
11	Draw over vaporizer - Isoflurane	9100-047
12	Draw over vaporizer - Halothane	9100-048
13	Breathing pressure gauge assembly	1500-005 + 2150-034 +1215-008
14	Reservoir bag	1300-070
15	Silicon breathing circuit	1300-082
16	Bellows	6650-022
17	Fenton balloon	6650-015
18	Balloon housing - clear part	2100-094
19	Patient Taper	2150-056
20	O ₂ single tube flowmeter	1525-009
21	Vaporizer back bar with 2 liter reservoir, air entrainment valve, over pressure valve, custom vaporizer connections	9000-062
22	Manual Ventilator Manifold Assembly (includes 2 unidirectional valves, 55 cmH ₂ O pressure relief valve and exhaust valve)	8150-162
23	Water trap and O seal	1300-068 + 1230-332
24	N ₂ O Cut Off Solenoid Assembly	8150-172
25	Exhaust non return valve and taper	2100-217

MONITORING

26	Control screen- includes touch screen and housing	NPN
27	Control screen circuit board	1400-009
28	Oxygen sensor	1300-075
29	Speaker	8150-179
30	Vital signs monitor	1419-221

ELECTRICAL

- 31 O₂ concentrator on/off switch
- 32 Fuses for O₂ concentrator and sockets - 5 amp x 5
- 33 Fuses for control screen - 500 Max 5
- 34 Main incoming fuse - 13 amp x 5
- 35 Mains inlet cable with UK style plug - 4 meters
- 36 Electrical sockets- 13 amp UK style
- 38 Mains isolator rocker switch
- 39 Power supply for control screen - 24V 15W
- 40 Automatic voltage switcher
- 41 Battery for control screen - YUASA NP1.2-12
- 12V 1.2 AH sealed lead acid battery

- 8150-174
- 1419-009
- 1419-010
- 1419-011
- NPN
- 1450-010
- 1450-009
- 1419-211
- 1450-015


- 1419-005



FRAME



- 42 Stainless steel work surface
- 43 Stainless steel drawer
- 44 Casters with brakes
- 45 Casters without brakes
- 46 Top monitor shelf


- 5150-093
- 5150-063
- 1300-092
- 1300-091
- 5150-067


ILLUSTRATED PARTS
REFER TO ITEM NUMBERS IN PARTS LIST

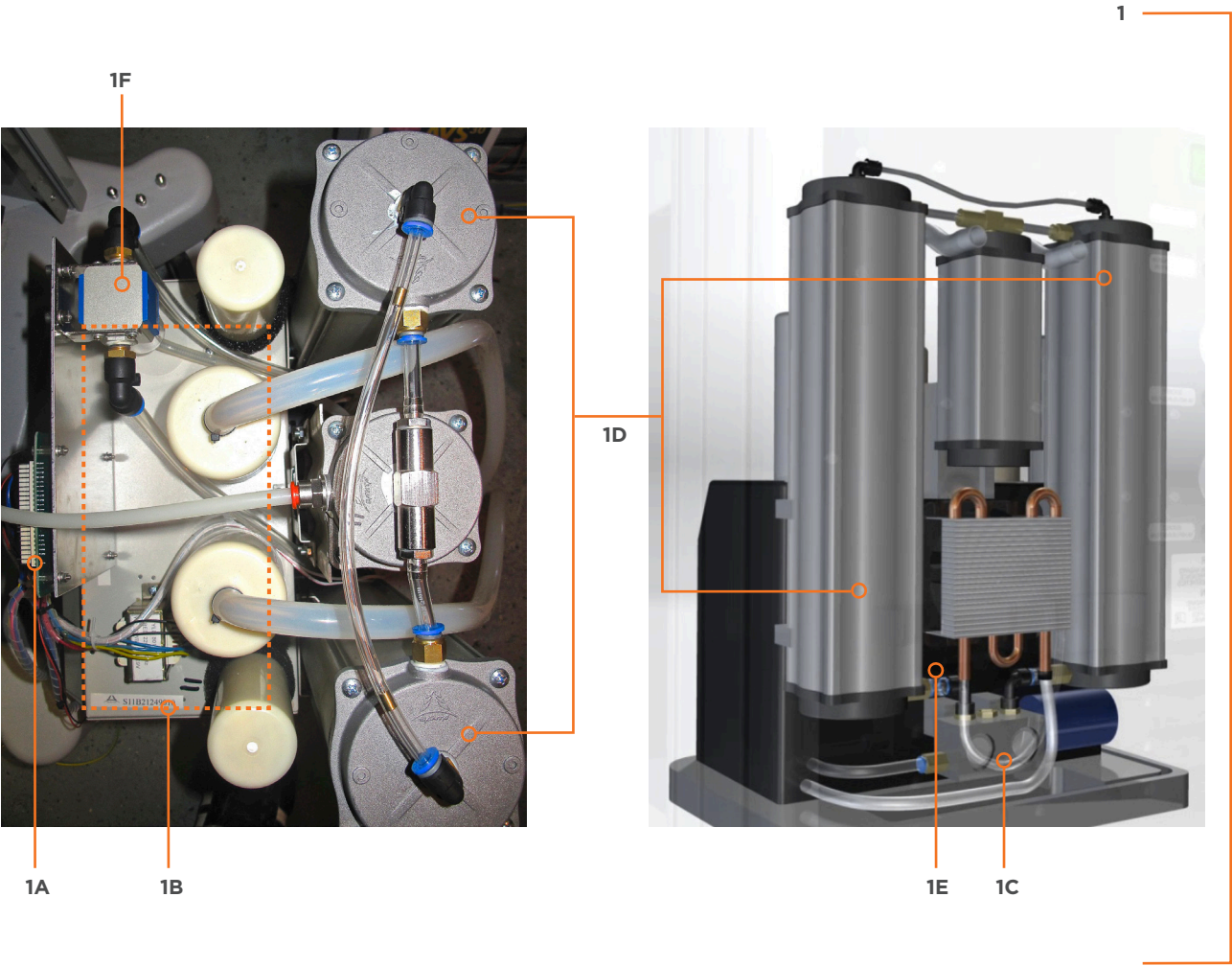

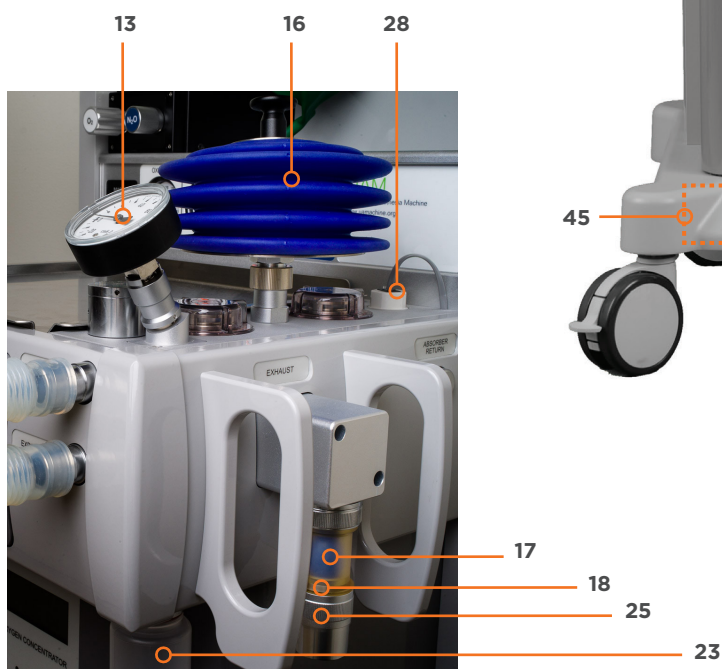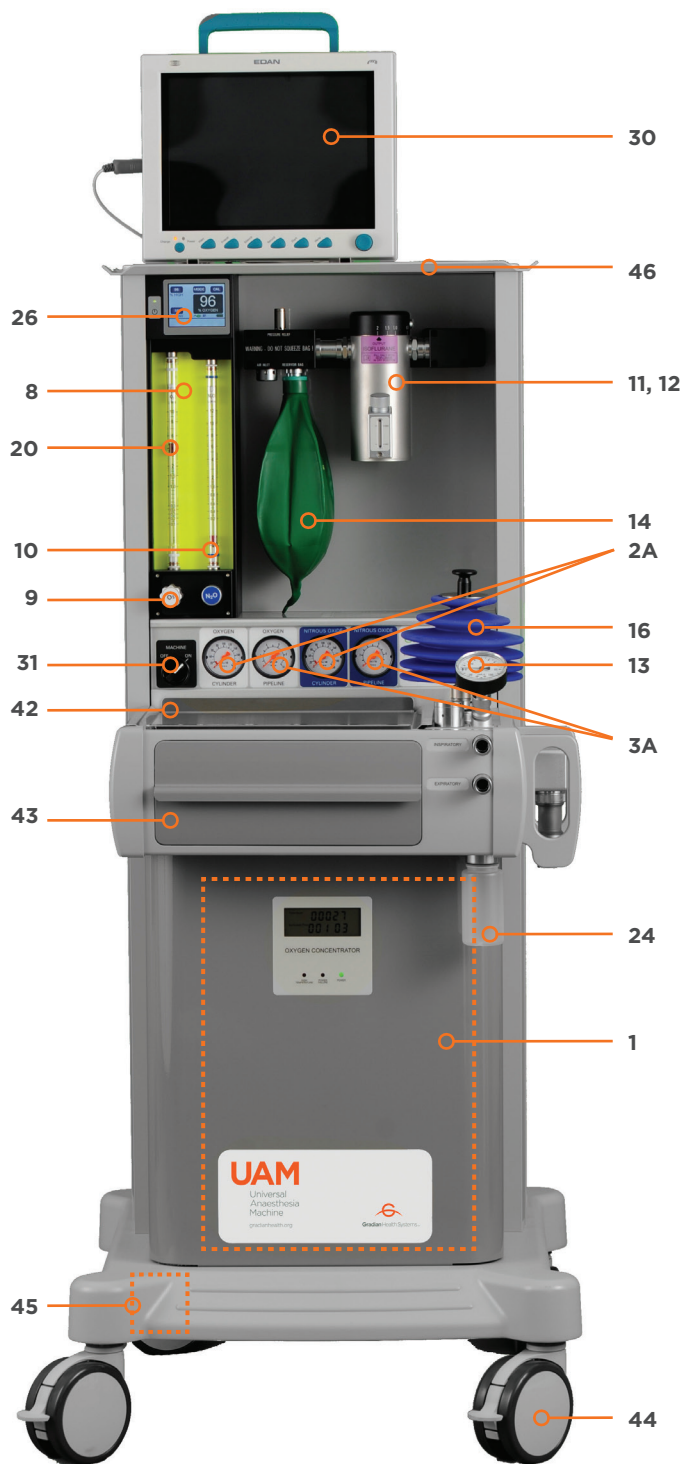


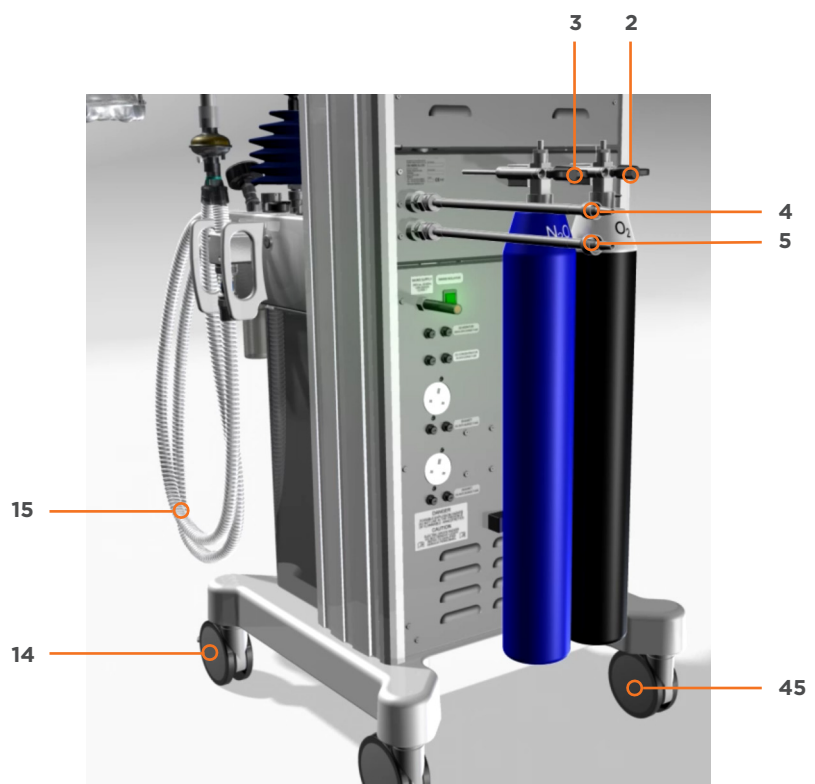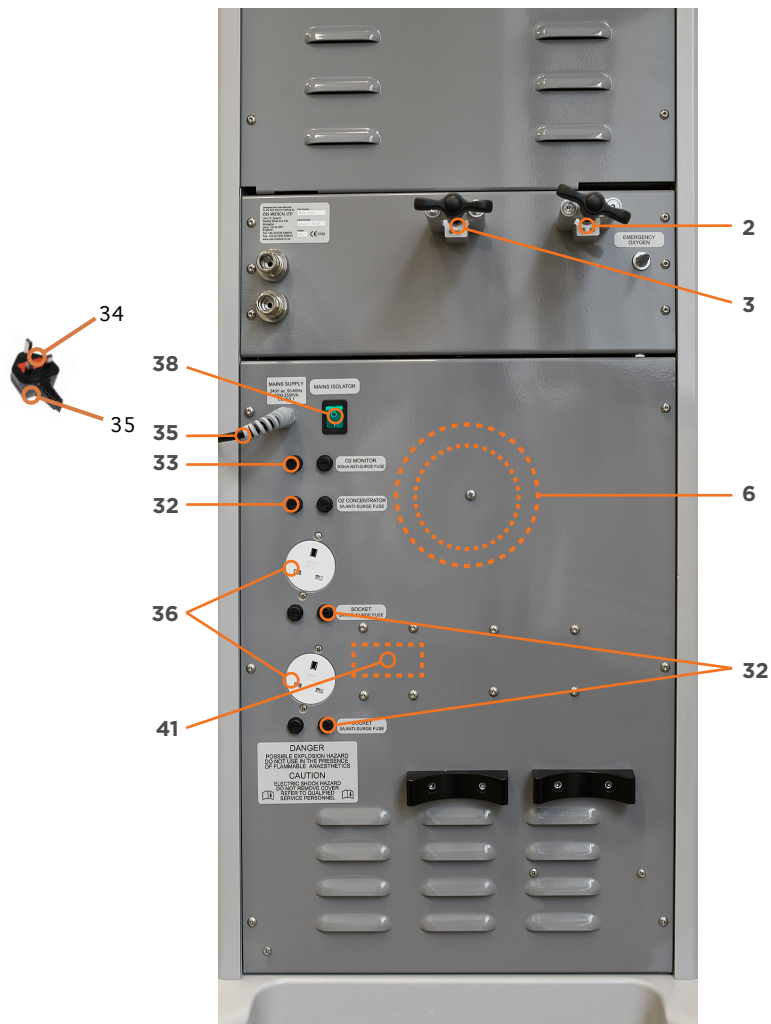

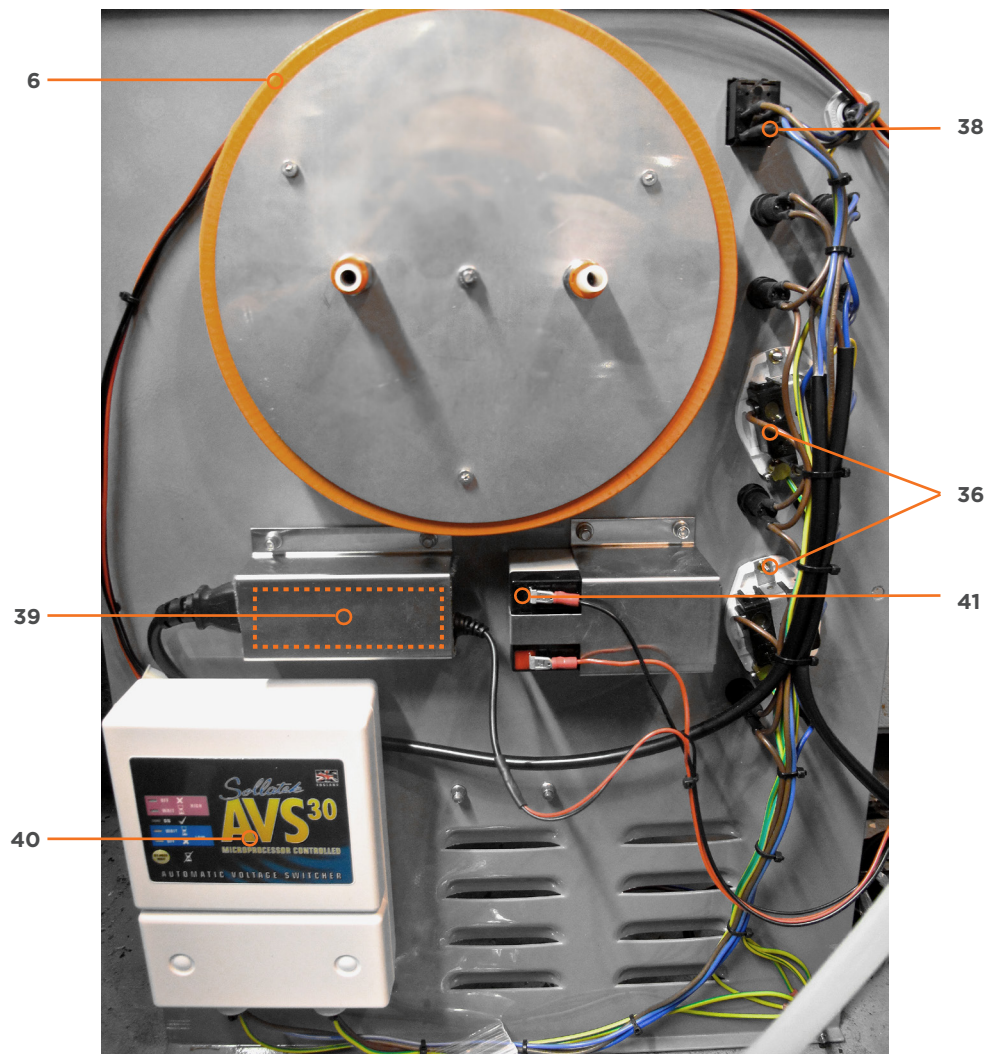


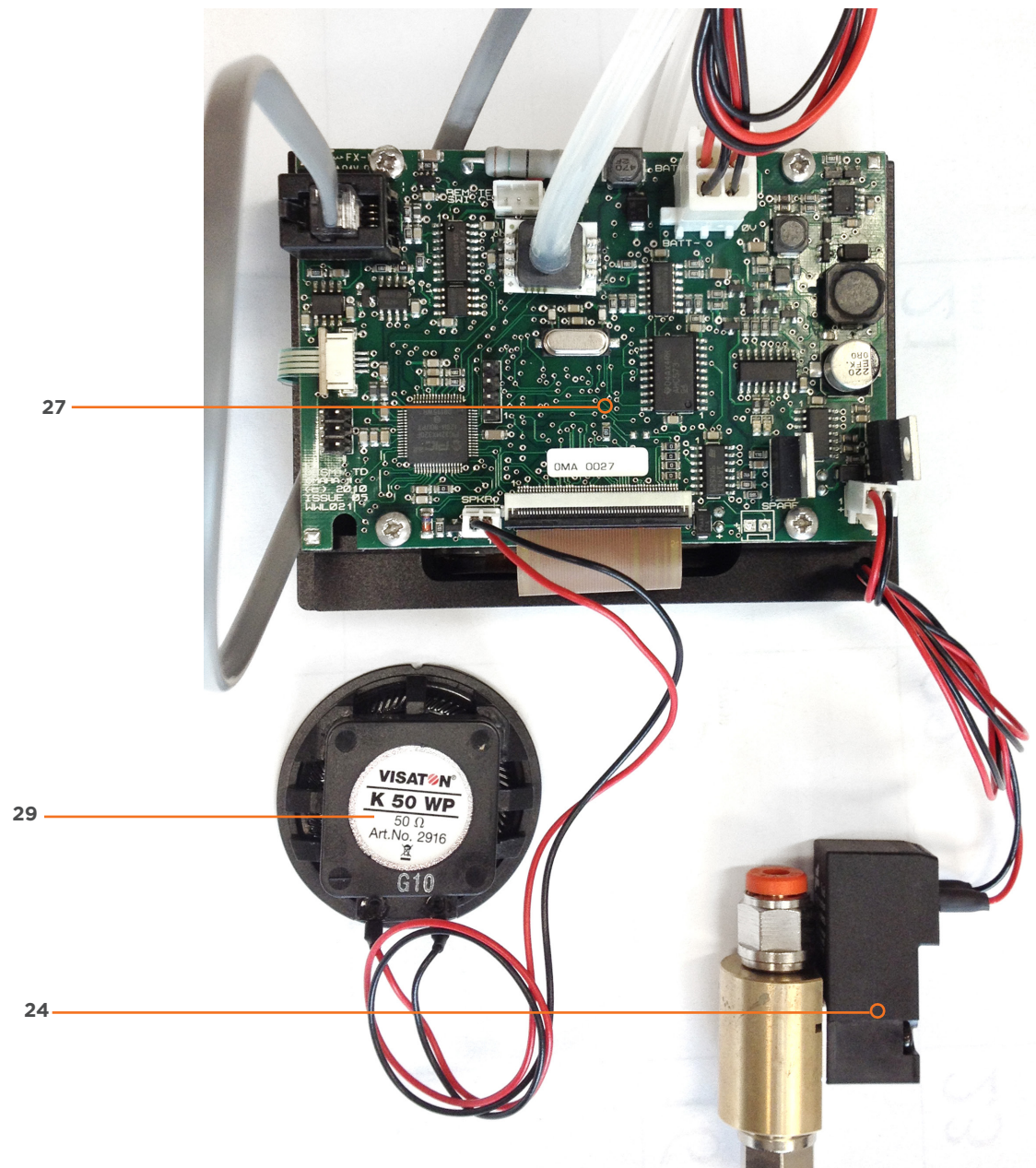


APPENDIX I
BREATHING CIRCUIT DIAGRAM

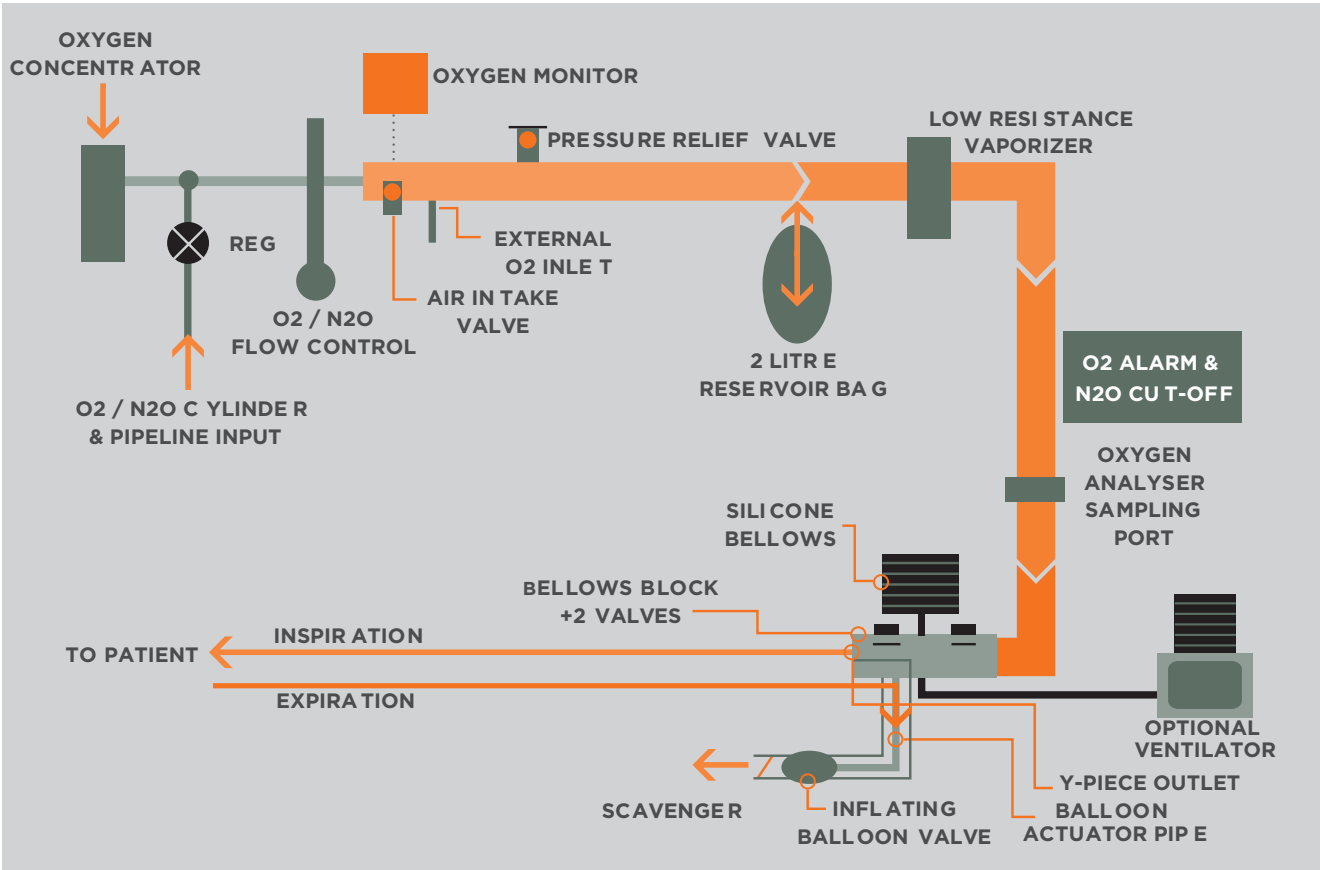

APPENDIX II
OXYGEN CONCENTRATOR FLOW DIAGRAM

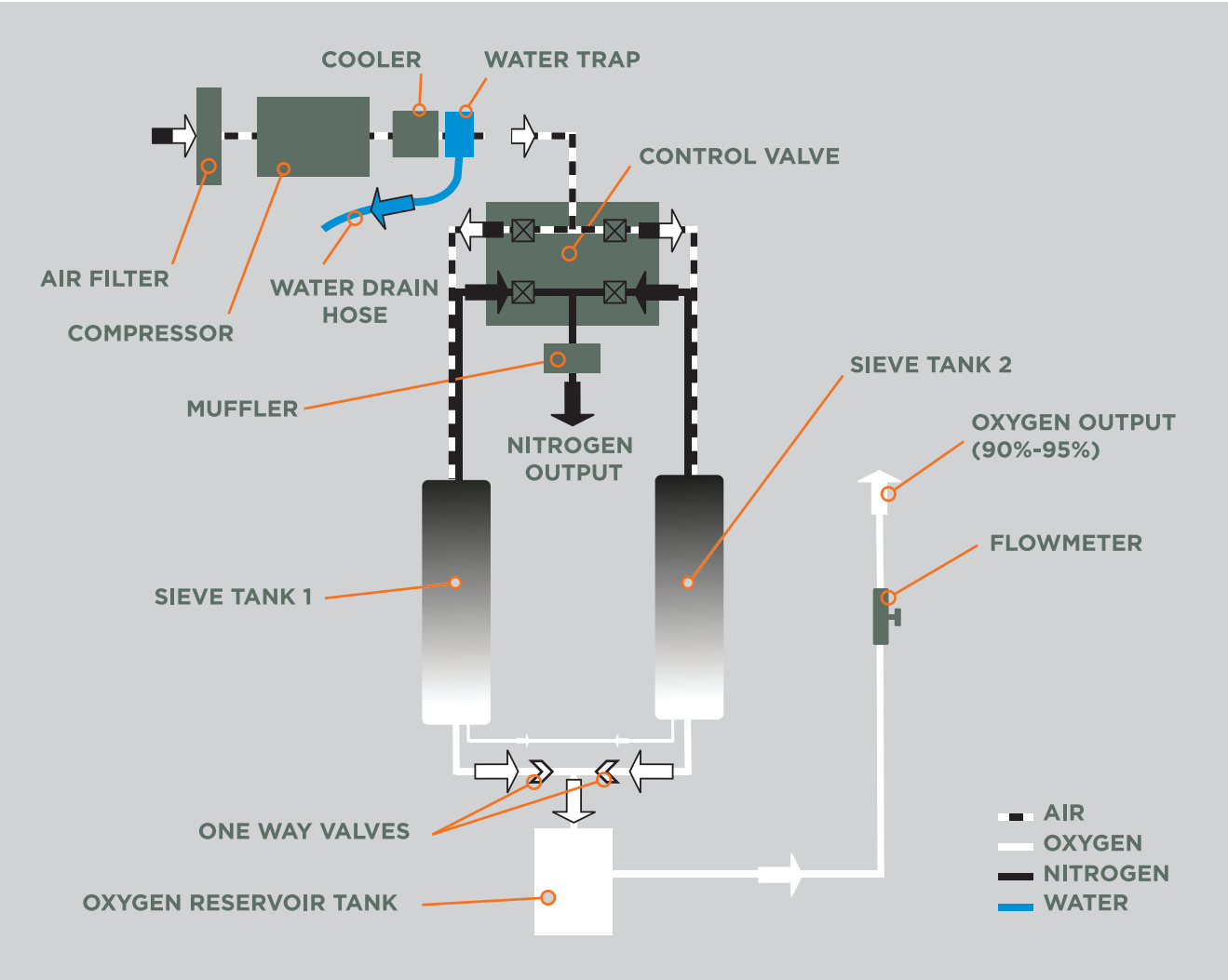

APPENDIX III
ELECTRICAL DIAGRAM

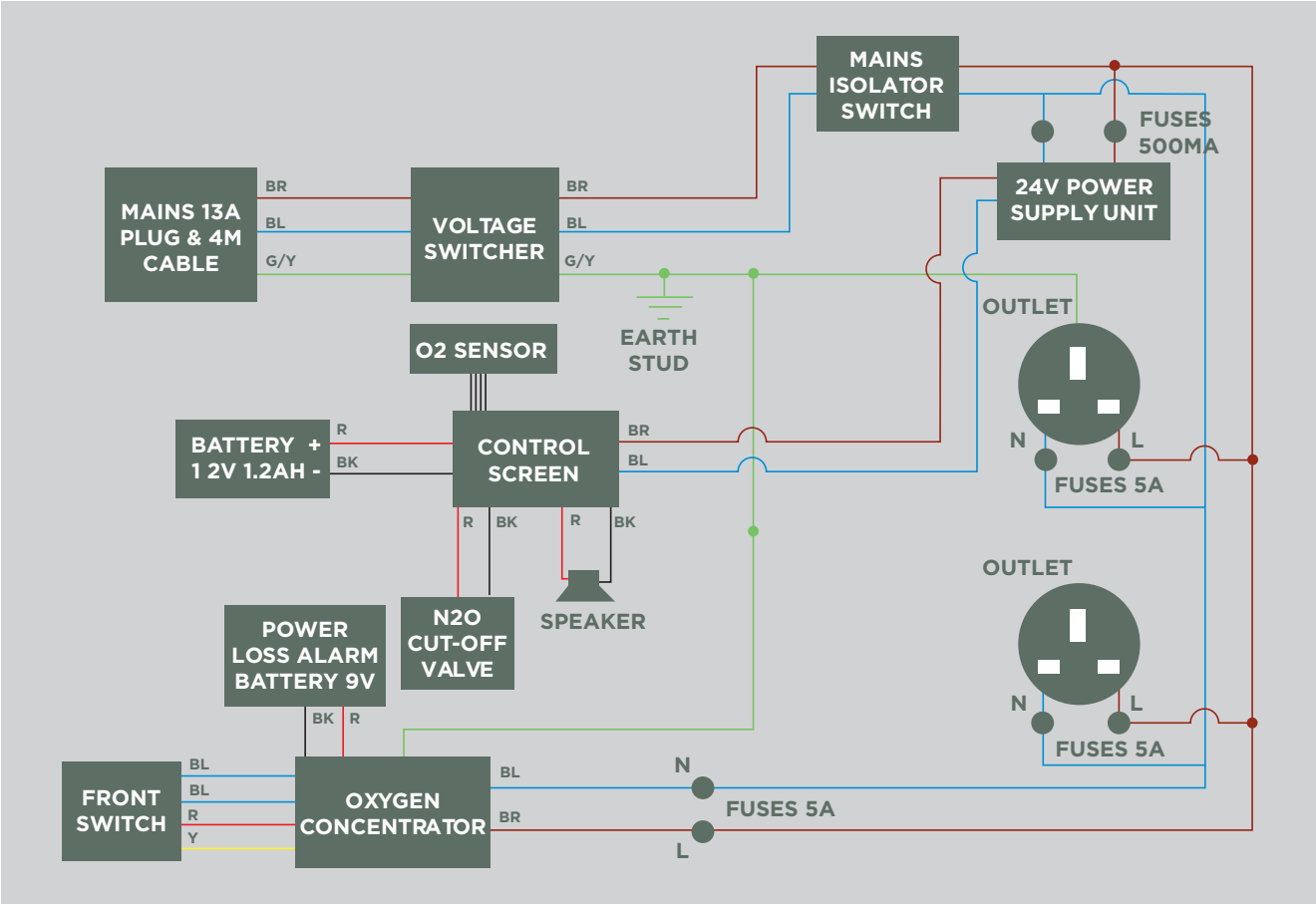

APPENDIX IV UAM PREVENTIVE MAINTENANCE AND FUNCTIONAL TESTS FORM

Carry these items when conducting preventive maintenance on a UAM: UAM Maintenance manual for instructions on PM and functional checks, new oxygen sensor, new Fenton balloon, new control screen backup 12 V battery, new 9 V battery, new air filter, new reservoir bag patient water trap. Only replace items if worn or defective.

MACHINE SERIAL # (located above the green mains isolator switch)		
CONCENTRATOR HOURS		
NEXT PREVENTIVE MAINTENANCE DUE DATE (6 months from today)		
CONTROL SCREEN SOFTWARE REVISION		
EXTERNAL GAS SOURCES Check and write down pressure reading. (If not available mark as N/A)		
OXYGEN CYLINDER		KPa
OXYGEN PIPELINE		KPa
NITROUS OXIDE CYLINDER		Kpa
NITROUS OXIDE PIPELINE		Kpa
SYSTEM	PASS OR FAIL	PROBLEM FOUND, WORK PERFORMED, PARTS USED
- All spares and accessories present & no modifications required		
- Frame, hardware, wheels		
- Mains cable and gas supply hoses		
- Charging indicator		
- Mains failure icon		
- Battery condition icon		
- Cylinder pin index yoke condition		
- Cylinder pressure gauges condition		
- Pipeline pressure gauges condition		
- Maximum O ₂ and N ₂ O flows (10 LPM)		
- Zero O ₂ and N ₂ O flows (0 LPM)		
- Concentrator air filter (clean first and (ONLY replace if cannot be cleaned)		
- Concentrator loss of power alarm battery (replace if no or low audio volume)		
OXYGEN CALIBRATION, SAFETY AND QUALITY	PASS OR FAIL	PROBLEM FOUND, WORK PERFORMED, PARTS USED
- Calibration with 100% O ₂ (use cylinder or pipeline oxygen)		
- Calibration with Air (use room air)		
- Upper and lower O ₂ alarms		
- N ₂ O cut-off		
- O ₂ sensor disconnect alarm		
- O ₂ concentrator output		
- O ₂ sensor (replace only if defective or will not calibrate properly)		
BREATHING SYSTEM	PASS OR FAIL	PROBLEM FOUND, WORK PERFORMED, PARTS USED
- Apnea alarm		
- Air intake valve function		
- 5 cmH ₂ O pressure relief valve function		
- Reservoir bag condition		
- Bellows assembly condition		
- Bellows expiratory pressure		
- Patient water trap condition		
- Fenton balloon operation (replace only if defective or has a hole)		
- 55 cmH ₂ O pressure relief valve		
VAPORIZER	PASS OR FAIL	PROBLEM FOUND, WORK PERFORMED, PARTS USED
- Vaporizer condition and indicator		
- Vaporizer performance		
OTHER	PASS OR FAIL	PROBLEM FOUND, WORK PERFORMED, PARTS USED
- Residual flow test		
- Perform ventilator checks if available (follow separate checklist)		
- Perform patient monitor checks if available (follow separate checklist)		
COMPLETE AND ADHERE SERVICE LABEL TO UAM		
Signature of Technician:		Hospital Personnel Signature:

APPENDIX V
UAM INSPECTION LABEL



Gradian Health Systems®

UAM INSPECTION LABEL

DATE: _____

BY: _____

RE-INSPECT BEFORE: _____

APPENDIX VI
PRESSURE CONVERSIONS

- Using the Table
1. Find the units you wish to convert FROM in the left hand Column.
 2. Find the units you wish to convert TO in the top row
 3. Insert the multiplier shown at the intersection into the following formula:

FROM units X MULTIPLIER = TO units
Example: 100 psi x 6.894757 = 689.475 kPa

DESIGNATION	PSI	KPA	CM OF H2O	BAR
PSI	1	6.894757	70.306927	.06894757
KPA	.1450377	1	10.19745	.01
CM OF H2O	.0142229	.0980634	1	.000980634
BAR	14.5038	100	1019.7466	1

COLOR IDENTIFICATION FOR MEDICAL GASES

GAS	ISO	USA
OXYGEN	White	Green
NITROUS OXIDE	Light Blue	Light Blue
MEDICAL AIR	Black and White	Yellow
SUCTION	Yellow	White

NOTES:



Gradian Health Systems®

Gradian Health Systems
915 Broadway, Suite 1001
New York, NY 10010, USA

gradianhealth.org

email: service@gradianhealth.org

tel: +1.212.537.0340

WhatsApp: +1.929.280.0210